

## Swift™ FX

### NASAL PILLOWS SYSTEM

#### Less mask. More you.

- Minimal design reinforces sense of freedom
- Unprecedented softness, simplicity and stability
- Flexible, fluid form follows facial contours and user movements
- Promotes early therapy acceptance
- Intuitive, virtually instant fitting
- Stable performance, even at high pressures

**Dual-wall nasal pillows**  
seal on contact for a gentle, comfortable seal, offer a wide fit range and disperse air gently to the nose

**Soft, narrow cushion with integrated flexible chamber**  
makes side sleeping comfortable and helps maintain seal during movement

**Ultra-smooth 360° rotating elbow**  
rotates with easy, fluid movement for flexibility in tube positions

**Quiet vents**  
diffuse and direct airflow quietly away from the user and bed partner

**Spring-flex lightweight tube**  
stretches and bends uniquely for increased stability without the use of a tube retainer

**Simple top buckle**  
has clearly marked notches for easy adjustment

**Soft silicone headgear**  
embraces the face gently and needs minimal adjustment—default setting fits most users


**Low-profile backstrap with Velcro fasteners**  
offers easy, fine adjustments

**Soft wraps**  
enhance comfort at the cheeks

**Size and name indicators**  
make it easy to identify the size as well as the correct orientation of the mask

**"Lock and key"**  
for easy attachment and detachment of side straps

**Flexible pillows with trampoline action**  
accommodate movement in multiple directions while maintaining seal

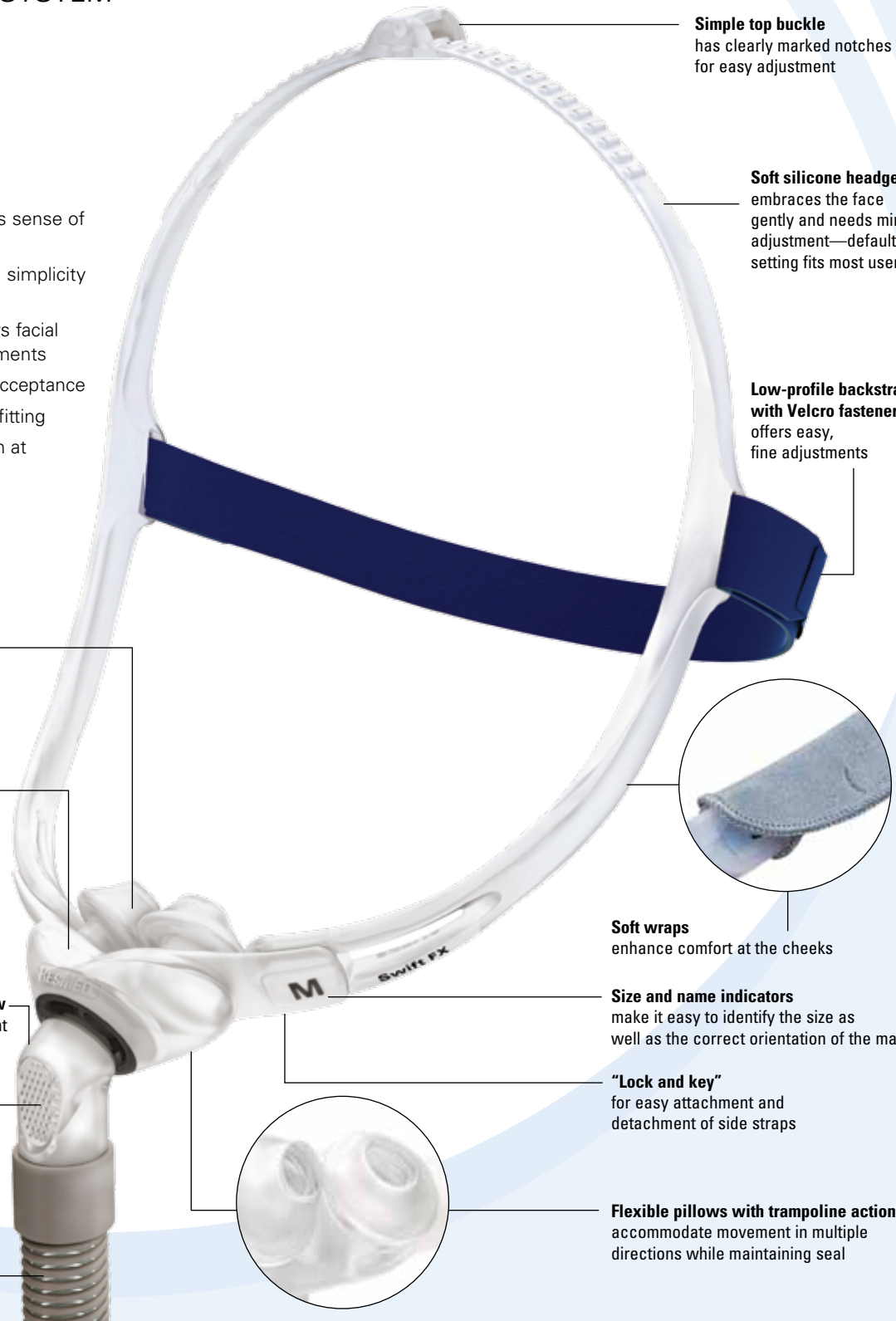



**MASK HIGHLIGHT**

**Comfort all around**

- With minimal facial contact and unobtrusive design, the Swift FX makes therapy less intimidating and easier to accept at an early stage.
- From the tip of the nasal pillows to the soft and secure sides, from the comforting soft wraps to the smooth backstrap ... the Swift FX is so light and gentle, patients will forget it's there.
- The snug headgear straps do the work of maintaining seal, so the fit feels loose and there's no need to over-tighten the headgear.

**USER BENEFITS**

**Soft comfort, stable performance**

- The soft design of the entire system enhances comfort.
- Dual-wall pillows and flexible cushion base maintain seal while adapting to movement.
- The flexible, lightweight tubing reduces pull on the mask and enhances stability.

**Patients may even forget it's there**

- Minimal design reinforces the sense of freedom—it's as much to do with how the user feels as with how the mask looks on their face.

**Choice of sleeping positions**

- Swift FX users can sleep on their side or their back, positioning the tube as desired.

**Quiet performance for quiet sleep**

- Innovative vent design delivers super-quiet mask performance—and gently disperses air away from patient and bed partner.

**PROFESSIONAL BENEFITS**

**Intuitive, almost instant fit**

- Saves valuable time for clinicians and sleep professionals. Once the right pillow size is selected, adjustment is minimal.

**Stable seal at high pressures**

- Uniquely designed, the system stays in place and maintains seal effectively, even at high treatment pressures.\*

**Easier therapy acceptance by patients**

- Designed for instant appeal, the Swift FX is a non-intimidating mask, facilitating early acceptance of therapy.

**Reduced inventory, wide (>95%) fit range\*\***

- Packed with three pillow sizes to suit most users, reducing inventory.
- An Extra Small size is also available.

**Easy product handling**

- With very few parts, the Swift FX is easy to clean and assemble.

**Convenient to recommend or prescribe**

- User guide and online support facilitate training for users and staff.



**PRODUCT CODES**  
**Swift FX Nasal Pillows System**

Americas 61500

\* Tested up to 20 cm H<sub>2</sub>O in clinical trials  
\*\* Based on ResMed's market research

# Swift™ FX

NASAL PILLOWS SYSTEM

## Components Card

### Complete system

61500 (includes Small (S), Medium (M) and Large (L) pillows\*)

### Frame system

(no headgear)

XS\* **61510**  
S **61511**  
M **61512**  
L **61513**

Headgear assembly  
(headgear, backstrap)  
**61529**

Short tube assembly  
(short tube, elbow and swivel)  
**61528 (1)**  
**61527 (10)**

Pillows  
**XS\* 61520**  
**S 61521**  
**M 61522**  
**L 61523**

Soft wraps  
**61530 (2)**

\* Extra Small (XS) pillows are available separately



## Swift™ FX

NASAL PILLOWS SYSTEM

### Fitting Guide



- Holding the pillows at your nose, pull the headgear over your head.
- Headgear should sit midway between the ears and eyes.



- Adjust the top strap first.
- Pull the strap through the buckle one notch at a time.
- Do not over-tighten.



- Adjust the backstrap to fine-tune the fitting if necessary.
- Do not over-tighten.
- The backstrap can be positioned over or under the hair.



- Swift FX pillows are designed to provide a natural fit.
- Minor adjustment can help position the pillows at the most comfortable angle in the nose.



- Connect the Swift FX short tube to the device air tubing.
- Swift FX is now fitted and ready for use.



- With the device turned on, check that the nasal pillows are properly inflated. To do this, pull the nasal pillows away from the face and cheek, make sure they are inflated and reposition them in the nose.

#### OPTIONAL SOFT WRAPS

- Wrap each of the soft wraps around either side of the headgear, folding the top fastener over the lower one.
- The soft wraps should be positioned on the cheeks, with the fasteners on the outside.



Note: For additional information, refer to the *Swift FX Nasal Pillows System User Guide*.