

RESMED



Click Here To Buy

Swift™ FX Nano

Small in size, big on comfort



The Swift FX Nano nasal mask combines the open, streamlined design of nasal pillows that patients love with the natural breathing comfort of a nasal mask.



Swift™ FX Nano
NASAL MASK

ResMed.com/SwiftFXNanoVideos

*Combined responses of those expressing a preference in an independent patient survey in the USA, UK, Germany, France & Australia, from September 2012. Visit www.resmed.com/maskbrand.



Swift FX Nano


Patients, especially those new to CPAP therapy, often look for the smallest, least intrusive mask system. The Swift FX Nano's built-in technology has enabled us to remove the rigid frames and cumbersome forehead support traditionally found in nasal masks without compromising performance. The result is a nasal mask that offers the best of both worlds: a minimalist design—similar to nasal pillows—that delivers reliable therapy for increased comfort and compliance.

Ideal for patients

Designed to have instant patient appeal, the Swift FX Nano is also easy to wear and clean, allows them to read, wear glasses or watch television, and helps them sleep comfortably through the night.

Ideal for your business

The Swift FX Nano's simple, compact design helps reduce patient anxiety at setup so they're more likely to stay compliant. And because it's incredibly easy to fit and maintain, you spend less time troubleshooting and answering patient callbacks.



Swift FX Nano for Her features a stylish headgear and small cushion size specifically for women

Quick fitting two-point headgear

with easy adjusting notches on the top strap that a patient can set once and forget

Optional soft wraps

provide added comfort and an extra defense against facial marks

Modified dual-wall cushion

with trampoline action provides a stable, comfortable seal to minimize mask leak

Ultra smooth ball-joint

elbow moves freely in virtually any angle and direction for a stable seal, reduced tube drag and greater sleeping freedom

Whisper-quiet venting

gently distributes air away from patient and bed partner

Spring-flex tubing and dual swivel

reduce pull on the mask to accommodate various sleeping positions

SoftEdge™ fabric backstrap

is soft and comfortable, with quick setting fasteners

Improved size and name indicators

clearly identify size and correct orientation of mask

"Lock and key" assembly

for easy attaching and detaching of headgear

Cushion tabs

provide total flexibility to compensate for additional movement during sleep



Swift FX Nano and Swift FX Nano for Her components

PRODUCT CODES

Complete System

Swift FX Nano Nasal Mask

Style	Part Number
Swift FX Nano (Standard) Standard cushion, Swift FX Nano headgear assembly, plus gray backstrap, soft wraps and frame assembly	62200
Swift FX Nano (Wide) Wide cushion, Swift FX Nano headgear assembly, plus gray backstrap, soft wraps and frame assembly	62251
Swift FX Nano for Her (Small) Small cushion, Swift FX Nano for Her headgear assembly, plus pink backstrap, soft wraps and frame assembly	62201

Reimbursement Codes

HCPCS	Code Descriptor
A7034	Mask interface (mask or cannula), 1 per 3 months
A7033	Pillow for use on nasal cannula type interface (replacement), 2 pairs per month
A7035	Headgear used with positive airway pressure device, 1 per 6 months

Headgear assembly


(headgear, backstrap and soft wraps)

62238 Gray
62248 Pink



Soft wraps

61530 Gray (2)
61544 Pink (2)

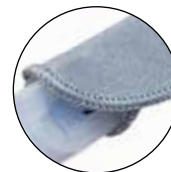

Cushions

62230 Standard
62281 Wide
62231 Small

Frame system

(no headgear)

62220 Standard
62271 Wide
62221 Small

Frame assembly (Gray)

(short tube, elbow and swivel)

62236 (1 pack)

Frame assembly (Pink)

(short tube, elbow and swivel)

62244 (1 pack)



Swift™ FX Nano
NASAL MASK

[ResMed.com/SwiftFXNanoVideos](https://www.resmed.com/SwiftFXNanoVideos)



Encourage your patients to join Wake Up to Sleep—ResMed's sleep apnea support community—at WakeUpToSleep.com/Join.

Swift FX Nano



- Ensure bottom of cushion opening aligns with the bottom of your nose, then place the mask over your nose.
- Then pull headgear over your head so it sits midway between ears and eyes.



- Adjust the top strap for positioning, one notch at a time, until the mask is comfortably over your nose. Do not overtighten.

TIP: If cushion is putting pressure on the bottom of nose, loosen top strap.



- The backstrap should sit just above the ears at the back of the head.
- Adjust both fasteners on the fabric backstrap at the same time.
- Ensure backstrap sits over your hair.



- Ensure mask is comfortably sealed.

TIP: When adjusting for leak, adjust backstrap



- Adjust the nasal cushion until it sits comfortably on your nose.
- Connect the free end of the short tube assembly to the air tubing of your device, then turn it on.



- If air is leaking, check for leak and adjust backstrap.
- Your mask is now ready to use.

TIP: For optimal performance, wear tube down



Optional Soft Wraps

- Wrap each of the soft wraps around either side of the headgear, folding the top fastener over the lower one.
- The soft wraps should be positioned on your cheeks, with the fasteners on the outside.

Note: For additional information, refer to the Swift FX Nano Nasal Mask User Guide.



Scan here to see fitting and care videos

www.resmed.com/SwiftFXNanoVideos



S9 Series
More Comfort.

Combine a Swift FX Nano mask with ResMed's S9 AutoSet™ therapy device for a comfortable, quiet and lifestyle-friendly system.

Swift™ FX Nano
NASAL MASK

ResMed.com/SwiftFXNanoVideos