

RESPIRONICS



 **RespShop**

**EverFlo™  
EverFlo™ Q**

USER MANUAL | MODE D'EMPLOI POUR LE PATIENT | MANUAL DEL PACIENTE | MANUAL DO PACIENTE

**RESPIRONICS**  
1001 Mario Ridge Lane  
Murtyville, PA 15668  
U.S.A.

EC REP

Respironics Deutschland  
Gemeinschaftsges. m.b.H.  
82031 Herrsching, Germany

CE  
0273













MS3274  
DS-41/08  
CH-INTL, PA-CA, ES-LAT, FR-08

# EverFlo / EverFlo Q User Manual

## Table of Contents

Symbol Key .....	2
Abbreviations .....	2
<b>Chapter 1: Introduction</b> .....	<b>3</b>
Intended Use .....	3
About Your EverFlo / EverFlo Q .....	3
Parts of Your Concentrator .....	3
Accessory Equipment and Replacement Parts .....	3
Warnings and Cautions .....	4
<b>Chapter 2: Operating Instructions</b> .....	<b>5</b>
<b>Chapter 3: Cleaning &amp; Maintenance</b> .....	<b>7</b>
Cleaning .....	7
Service .....	7
How to Contact Respironics .....	7
<b>Chapter 4: Alarms and Troubleshooting</b> .....	<b>8</b>
Troubleshooting Guide .....	9
<b>Chapter 5: Specifications</b> .....	<b>10</b>
<b>Appendix A: EMC Information</b> .....	<b>12</b>
<b>Limited Warranty</b> .....	<b>14</b>

## Symbol Key

	Follow Instructions for Use		AC Power
	Type of Applied Part	<b>REF</b>	Model Number
	Class II equipment	<b>SN</b>	Serial Number
	No smoking		On (Power)
	No oil or grease	○	Off (Power)
	Do not disassemble	<b>IPX1</b>	Drip proof equipment
	General Alarm		Action Required, Check System Indicator
	Canadian/US Certification		Compliant with the Waste Electrical and Electronic Equipment/Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment (WEEE/RoHS) recycling directives
	European Declaration of Conformity		

## Abbreviations

LED	Light Emitting Diode
LPM	Liters per Minute
OPI	Oxygen Percentage Indicator

ResFlo is a trademark of Respiration, Inc. and its affiliates.  
 © 2008 Respiration, Inc. and its affiliates. All rights reserved.

## Chapter 1: Introduction

Your health care professional has determined that supplemental oxygen is of benefit to you and has prescribed an oxygen concentration set at a specific flow setting to meet your needs. **DO NOT** change the flow settings unless your health care professional tells you to do so. **Please** read and understand this entire manual before using the device.



### Intended Use

The EverFlo / EverFlo Q Oxygen Concentrator is intended to provide supplemental oxygen to persons requiring oxygen therapy. The device is not intended to be life supporting or life sustaining.

### About Your EverFlo / EverFlo Q

The device produces concentrated oxygen from room air for delivery to a patient requiring low flow oxygen therapy. The oxygen from the air is concentrated using a molecular sieve and a pressure swing adsorption process. Your home care provider will show you how to operate the concentrator and will be available to answer any questions. If you have additional questions or problems, contact your home care provider.

### Parts of Your Concentrator

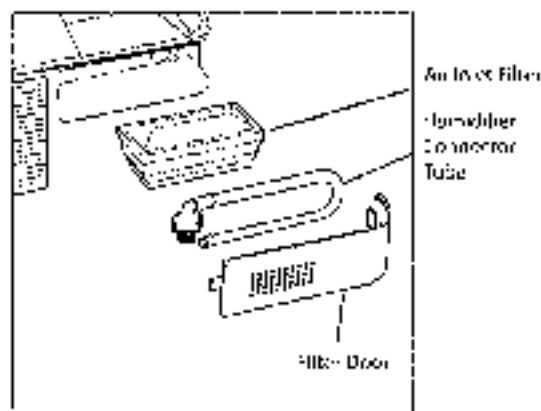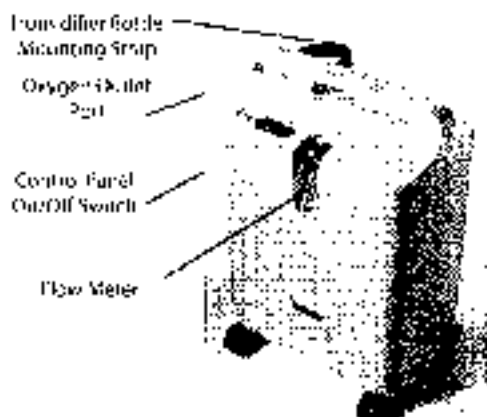

### Accessory Equipment and Replacement Parts

Contact your home care provider if you have questions about this equipment. Use only the following Respiration accessories and replacement parts with this device:

- Air Inlet Filter
- Humidifier Connector Tube

## Warnings and Cautions

### Warnings

*A warning represents the possibility of harm to the operator or patient.*

- For proper operation, your concentration requires unobstructed ventilation. The ventilation ports are located at the rear base of the device and at the side inlet filter. Keep the device at least 15 to 30 cm away from walls, furniture, and especially curtains that could impede adequate airflow to the device. Do not place the concentrator in a small closed space (such as a closet).
- Do not remove the covers of this device. Removing covers is referred to as authorized maintenance. Respiratory home care involves.
- In the event of an equipment alert, and if you are experiencing any signs of discomfort, consult your home care provider and/or your healthcare professional immediately.
- Oxygen generated by this concentrator is supplemental and should not be considered life supporting or life sustaining. In certain circumstances oxygen tanks can be hazardous. Any user should seek medical advice prior to using this device.
- Where the prescribing health care professional has determined that an interruption in the supply of oxygen, for any reason, may have serious consequences to the user, an alternate source of oxygen should be available for immediate use.
- Oxygen vigorously accelerates combustion and should be kept away from heat or open flames. Not suitable for use in the presence of flammable anesthetic mixture with an oxidizing agent or nitrous oxide.
- Do not smoke, allow others to smoke or have open flames near the concentrator when it is in use.
- Do not use oil or grease on the concentrator or its components as these substances when combined with oxygen, can greatly increase the potential for a fire hazard and personal injury.
- Do not use the oxygen concentrator if either the plug or power cord is damaged. Do not use extension cords or essential adapters.
- Do not attempt to clean the concentrator while it is plugged into an electrical outlet.
- Device operation above or outside of the voltage, LPM, temperature, humidity and/or altitude values specified may decrease oxygen concentration levels.
- Your home care provider is responsible for performing appropriate preventive maintenance at the intervals recommended by the device manufacturer.

### Cautions

*A caution represents the possibility of damage to the equipment.*

- Do not place liquids on or near the device.
- If liquid is spilled on the device, turn the power off and unplugging from electrical outlet before attempting to clean up spill. Call your home care provider if device does not continue to work properly.


## Chapter 2: Operating Instructions

**Warning:** Do not use extension cords or electrical adapters.

1. Select a location that allows the concentrator to draw in room air without being restricted. Make sure that the device is at least 15 to 30 cm away from walls, furniture, and especially curtains that could impede adequate air flow to the device. Do not place the device near any heat source.
2. After reading this entire manual, plug the power cord into an electrical outlet.
3. Do either Step A or Step B below.
  - A. If you are not using a humidifier, connect your nasal cannula to the Oxygen Outlet Port, as shown in the next illustration on the right.
  - B. If you are using a humidifier, follow the steps below:
    1. Open the filter door on the back of the device as shown.
    2. Remove the humidifier connector tube from the back of the filter door and replace the filter door, as shown.
    3. Fill your humidifier bottle according to the manufacturer's instructions.
    4. Place the filter humidifier on the top of the EverFlo / EverFlo Q device inside the velcro strap, as shown in the illustration on the right.
    5. Tighten the velcro strap around the bottle and secure it so it is held firmly in place.
    6. Connect the humidifier connector tube (that you received from the filter door) to the Oxygen Outlet Port (as shown in Step 3-A above).
    7. Connect the other end of the humidifier connector tube to the top of the humidifier with the elbow in the wrong facing the Front, as shown here.
    8. Connect your cannula to the humidifier bottle according to the humidifier bottle manufacturer's specifications.




Step 3-B1

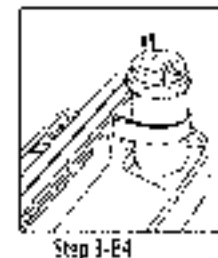


## EverFlo / EverFlo Q User Manual

4. Press the power switch to the On (I) position. Initially, all the LEDs will illuminate and the audible alert will beep for a few seconds. After that time, only the green LED should remain lit. You can begin breathing from the device immediately even though it typically takes 10 minutes to reach oxygen purity specifications.

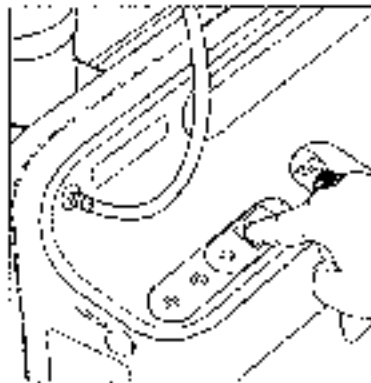

5. Adjust the flow to the prescribed setting by turning the knob on the top of the flow meter until the ball is centered on the line marking the specific flow rate.

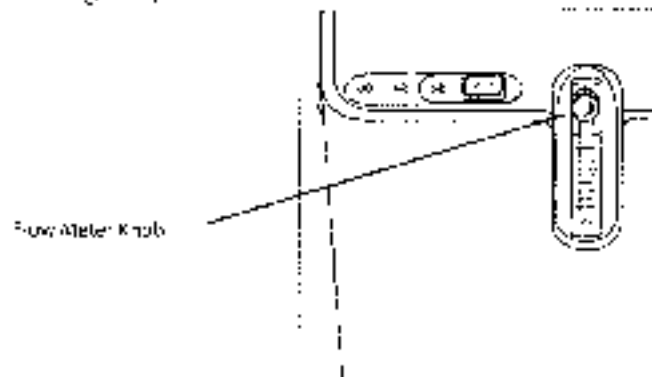

6. Be sure oxygen is flowing through the cannula. If it is not, refer to the User's Adjusting Guide in this manual.
7. Put on the cannula as directed by your home care provider.
8. When you are not using the oxygen concentrator, press the power switch to the Off (O) position.

## Chapter 3: Cleaning & Maintenance

*Warning: It is important to unplug the device before you perform any cleaning.*

*Caution: Excess moisture may impair the proper operation of the device.*

### Cleaning

Periodically, use a damp cloth to wipe down the exterior case of the EverFlo / EverFlo Q device. If you use medical disinfectants, be sure to follow the manufacturer's instructions.

If you are using a humidifier, clean your device according to your home care provider's or manufacturer's instructions.

### Service

The EverFlo / EverFlo Q Oxygen Concentrator contains no user-serviceable parts.

*Warning: Do not remove the cover of this device. Servicing must be referred to an authorized and trained Respironics home care provider.*

### How to Contact Respironics

To have your device serviced, contact your home care provider. If you need to contact Respironics directly, call the Respironics Customer Service department at 1-724-587-4000 or Respironics Deutschland at +49 8152 93060. You can also use the following addresses:

Respironics  
1001 Marty Ridge Lane  
Murrysville, PA 15668  
USA

Respironics Deutschland  
Gewerbestraße 17  
82211 Herrsching  
Germany


Visit the EverFlo website at [www.everflo.respironics.com](http://www.everflo.respironics.com)



## Chapter 4: Alarms and Troubleshooting

### Alarm and Indicators

The device has an audible alarm and three LED indicators, as shown below.

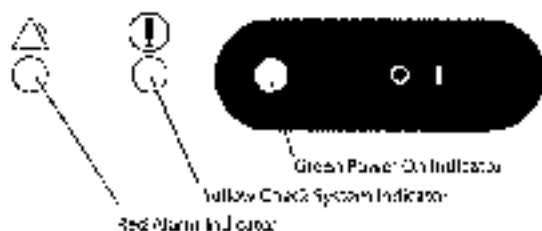

Audible Alarm / Colored LED	Possible Cause	Your Action
All 3 LEDs illuminate continuously and the Audible Alarm is sounding continuously.	The device has detected a system malfunction.	Immediately turn off the device, connect to a back up oxygen source, and call your home care provider.
The Audible Alarm is sounding continuously. None of the LEDs are illuminated.	The device is turned off but is not operating. Often this indicates that the device is not plugged in or there is a power failure.	Check the power outlet and verify that the device is plugged in. If the problem continues, connect to a back up oxygen source and call your home care provider.
Red LED illuminates continuously and the Audible Alarm is sounding continuously.	The device has detected a system malfunction.	Immediately turn off the device and wait 3 minutes. Restart the device. If the condition persists turn the unit off, connect to a back up oxygen source, and call your home care provider.
Yellow LED illuminates continuously. The Red LED is blinking and the Audible Alarm is beeping periodically.	The device has detected an impeded oxygen flow condition.	Follow the troubleshooting guide on the next page. Connect to a back up oxygen source and call your home care provider if your troubleshooting efforts fail to end this alert condition.
Red LED illuminates continuously. The Red LED is off and the Audible Alarm is silent.	The device has detected a low oxygen condition (OP) units only.	Continue using the unit but call your home care provider about this condition.

## Troubleshooting Guide

Problem	Why it Happened	What to Do
Yellow LED is blinking. The Red LED is off and the Audible Alarm is beeping periodically.	The device has detected a high-oxygen flow condition.	Turn the flow rate down to your prescribed level. Wait at least 2 minutes. If the condition persists turn the unit off, connect to a back up oxygen source, and call your home care provider.
Green LED illuminates continuously. The other LEDs are off and the Audible Alarm is silent.	The device is turned on and working properly.	Take no action.
The device is not working when it is turned on.  (The Audible Alarm is sounding continuously. All LEDs are off.)	The power cord plug is not properly inserted into the electrical outlet.	Make sure the device is properly plugged in to the electrical outlet.
	The unit is not receiving power from the electrical outlet.	Check your household outlet fuse or circuit.
	Internal part failure.	Connect to a back up oxygen source and contact your home care provider.
The device is not working when it is turned on.  (The Audible Alarm is sounding continuously and all 3 LEDs are illuminated.)	Internal part failure.	Connect to a back up oxygen source and contact your home care provider.
Impeded oxygen flow indication is activated.  (The Yellow LED illuminates continuously, the Red LED is blinking and the Audible Alarm is beeping periodically.)	The airflow to the device is impeded or blocked.	Remove any items that appear to be blocking the airflow into the device.
	The flow meter knob is completely closed.	Turn the flow meter knob counterclockwise to center the ball on the prescribed LPM flow.
	The oxygen tubing is kinked and blocking the delivery of oxygen.	Check to see that the tubing is not kinked or blocked. Replace if necessary.
Limited oxygen flow to the user without any fault indicator.  (All LEDs are off and the Audible Alarm is silent.)	The oxygen tubing or cannula is faulty.	Inspect and replace the items if necessary.
	There is a poor connection to a device accessory.	Ensure that all connections are free from leaks.

## Chapter 5: Specifications

### Environmental

	Operating	Transport & Storage
<b>Temperature</b>	13 to 32°C	-11 to 71° C
<b>Relative Humidity</b>	15 to 95%, non-condensing	15 to 95%, non-condensing
<b>Altitude</b>	0 to 2286m	N/A

### Physical

<b>Dimensions</b>	58 cm x 33 cm x 24 cm
<b>Weight</b>	14 to 15 kg

### Standards Compliance

This device is designed to conform in the following standards:

- IEC 60601-1 Medical Electrical Equipment, Part 1: General Requirement for Safety
- IEC 60601-1-2 2nd edition, Medical Electrical Equipment, Part 1-2: General Requirement for Safety - Collateral standard, Electromagnetic Compatibility - Requirements and Tests
- ISO 3359 Oxygen Concentrators for Medical Use - Safety Requirements

### Electrical, AC Power Consumption

<b>Models 1020000, 1020001 1020002, 1020003 1020014, 1020015</b>	120 VAC ±10%, 350 W, 60 Hz
<b>1020004, 1020005</b>	230 VAC ±10%, 120 W, 60 Hz
<b>1020006, 1020007, 1020008 1020009, 1020010 1020011, 1020012 1020016, 1020017</b>	230 VAC ±10%, <300 W, 50 Hz
<b>1020013</b>	230 VAC ±10%, <300 W, 60 Hz

## Oxygen

<b>Oxygen Concentration*</b> (All Models except 1020007, 1020008)	50-96% from 0.5 to 51.2% <sup>†</sup>
<b>Models 1020007, 1020008</b>	87-96% from 0.5 to 51.2% <sup>†</sup>

- \* Device operation above or outside of the voltage, LPM, temperature, humidity or flow rate values specified may decrease oxygen concentration levels.

## Sound Level

<b>Models 1020000, 1020001</b>	45 dBA typical
<b>1020002, 1020003</b>	
<b>1020013</b>	
<b>1020004, 1020005</b>	43 dBA typical
<b>1020006, 1020008</b>	
<b>1020009, 1020010</b>	
<b>1020011, 1020012</b>	
<b>1020016, 1020017</b>	
<b>1020007, 1020014</b>	<40 dBA typical
<b>1020015</b>	

## Classification

The EverFlo / EverFlo Q Oxygen Concentrator is classified as:

- IEC Class II Equipment
- Type BF Applied Part
- LXI Drop Product
- Not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- Continuous Operation

## Disposal

Dispose of the device in accordance with local regulations.

## WEEE/RoHS Recycling Directives

If you are subject to the WEEE/RoHS recycling directives, refer to [www.respirics.com](http://www.respirics.com) for the program for recycling this product.

## Appendix A: EMC Information

**GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC EMISSIONS:** This device is intended for use in the electromagnetic environment specified below. The use of this device should make sure it is used in such an environment.


Emissions Test	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT – Guidance
RF emissions CISPR 11	Group 1	The device uses AC energy only for its internal function. Therefore, its RF emissions are very low and do not likely to cause any interference in nearby electronic equipment. This device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Class A	

**GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY:** This device is intended for use in the electromagnetic environment specified below. The use of this device should make sure it is used in such an environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT – Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	+6 kV contact -8 kV air	±6 kV contact ±8 kV air	Fluorescent light fixtures, overhead or portable life, floors are covered with synthetic material, the distance between floors should be at least 10cm.
Electrical Fast Transient Burst IEC 61000-4-4	+2 kV for power supply lines ±1 kV for input/output lines	±2 kV for supply lines ±1 kV for input/output lines	Main power quality should be that of a typical home or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV for common mode	Main power quality should be that of a typical home or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	+95% U <sub>n</sub> -95% dip in U <sub>n</sub> for 0.5 cycle 0% U <sub>n</sub> 100% dip in U <sub>n</sub> for 5 cycles 70% U <sub>n</sub> 130% dip in U <sub>n</sub> for 25 cycles rate U <sub>n</sub> 10-95% dip in U <sub>n</sub> for 5 sec	95% U <sub>n</sub> 100% dip in U <sub>n</sub> for 0.5 cycle 80% U <sub>n</sub> 60% dip in U <sub>n</sub> for 5 cycles 90% U <sub>n</sub> 130% dip in U <sub>n</sub> for 25 cycles rate U <sub>n</sub> 0-95% dip in U <sub>n</sub> for 5 sec	Main power quality should be that of a typical home or hospital environment. If the use of the device requires continuous operation during power mains interruptions, it is recommended that the device be powered from an independent power source or a battery.
Power frequency 50/60 Hz magnetic field IEC 61000-4-3	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical hospital or home environment.

Note: U<sub>n</sub> is the a.c. mains voltage prior to application of the test level.

**GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY:** This device is intended for use in a high electromagnetic immunity environment. The use of this device should only be in its intended and approved environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Conducted RF IEC 61024-4-6	3Vrms 30 kHz to 80 MHz	3Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 60601-4-3	1V/m 80 MHz to 2.5 GHz	1V/m	<p>Recommended separation distance:</p> $d = 1.2 \sqrt{P}$ $d = 1.2 \sqrt{P} \quad 30 \text{ MHz to } 800 \text{ MHz}$ $d = 2.2 \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ <p>Where P is the maximum output power rating of the transmitter (in Watts) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey<sup>a</sup>, should be no less than the compliance level in each frequency range<sup>b</sup>.</p> <p>Interference may occur if the vicinity of equipment marked with the following symbol  is used.</p>

Note 1: At 30 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a: Field strengths from fixed transmitters, such as base stations for mobile telephony, radio broadcasting, mobile satellite systems, FM and TV broadcast stations and IT systems, cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the estimated field strength at the location in which the device is used exceeds the applicable RF compliance level above, the device might be expected to be in non-compliance. If the normal performance is observed, additional measures may be necessary, such as re-orienting or re-locating the device.

b: Over the frequency range 150 kHz to 80 kHz, the field strength should be less than 0.1 V/m.

**RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THIS DEVICE:** The device is intended for use in a high electromagnetic environment in which cellular RF fields are common and controlled. The use of the use of this device can help prevent the absorption of interference by the user. The following table provides portable and mobile RF communications equipment (transmitters) and this device is recommended to be used according to the maximum output power of the communications equipment.

Rated Maximum Power Output of Transmitter (W)	SEPARATION DISTANCE According to Frequency of Transmitter (m)		
	150 MHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.2 \sqrt{P}$
0.01	0.12	0.12	0.22
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance (d in meters (m)) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in Watts (W) according to the transmitter manufacturer.

Note 1: At 30 kHz and 800 kHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

## Limited Warranty

Respironics, Inc. warrants that the system shall be free from defects of workmanship and materials and will perform in accordance with the product specifications for a period of three (3) years from the date of sale by Respironics, Inc. to the dealer. Respironics warrants that the EverFlo/ EverFlo Q units serviced by Respironics, or an authorized service center, will be free from defects in serviced materials for a period of 90 days and free from defects in workmanship for a period of 90 days from the time of service. Respironics accessories are warranted to be free of defects in materials and workmanship for a period of 90 days from the time of purchase. If the product fails to perform in accordance with the product specifications, Respironics, Inc. will repair or replace – at its option – the defective material or part. Respironics, Inc. will pay customary freight charges from Respironics, Inc. to the dealer location only. This warranty does not cover damage caused by accident, misuse, abuse, alteration, and other defects not related to material or workmanship.

Respironics, Inc. disclaims all liability for consequential loss of profits, overhead, or consequential damages which may be claimed to arise from any sale or use of this product. Some states do not allow the exclusion or limitation of incidental or consequential damages, so this above limitation or exclusion may not apply to you.

This warranty is given in lieu of all other express or implied warranties including the implied warranties of merchantability and fitness for a particular purpose. In addition, in no event shall Respironics be liable for lost profits, loss of good will, or incidental or consequential damages, even if Respironics has been advised of the possibility of the same. Some states or provinces do not allow the exclusion or limitation of implied warranties or the disclaimer of incidental and consequential damages. Accordingly, the laws of your state or province may give you additional protections.

To exercise your rights under this warranty, contact your local authorized Respironics, Inc. dealer or contact Respironics, Inc. at:

Respironics  
1031 Murry Ridge Lane  
Murrysville, PA 15668  
USA  
1-724-387-4000

Respironics Deutschland  
Gewerbestraße 17  
82211 Herrsching  
Germany  
+49 8152 93060