

# It's as easy as 1, 2, 3...

## 1. Fill

Simply lift the green lid and pour in the medication. Close the lid and make sure your child is sitting up straight, ready for therapy.



## 2. Treat

Press the on button and encourage your child to breathe in and out through the mouthpiece.\*\*\* InnoSpire Go will indicate therapy completion with audible and visible signals – before turning itself off.




## 3. Clean

After each treatment remove the mouthpiece and rinse it in warm water. At the end of each day wash it in warm soapy water and once a week disinfect it in boiling water.



Find out more

[www.philips.com/innospirego](http://www.philips.com/innospirego)

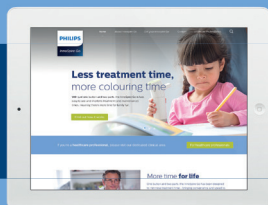

\*\* Using 2.5ml salbutamol

\*\*\* A mask is included for children that find it difficult to use the mouthpiece.



© 2016 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Respireonics Inc. and Respireonics Respiratory Drug Delivery (UK) Ltd are Philips companies.

Respireonics Respiratory Drug Delivery (UK) Ltd  
Chichester Business Park  
City Fields Way  
Tangmere  
Chichester  
PO20 2FT UK

KW 09/09/2016 PN452287  
PN 1132089 MCI 4107771  
[www.philips.com/respironics](http://www.philips.com/respironics)

**PHILIPS**

*InnoSpire Go*

Portable  
mesh nebulizer



More  
time  
for life

# More time for the yummy things in life

Fast, effective and easy to use, the InnoSpire Go nebulizer is specifically designed to shorten asthma treatment time so your child can get back to doing something fun.

## Reduces treatment time by 25%\*



### Treatment in 4 minutes\*\*

An efficient asthma therapy device that indicates therapy completion using audible and visible signals - before automatically turning off.



### Virtually silent

A small, portable and discreet nebulizer that can be used anywhere.



### Built-in long life battery


The rechargeable battery provides up to 30 treatments between each charge.



### Simple - only two parts

The easily detachable mouthpiece makes the InnoSpire Go simple to use, clean and maintain.\*\*\*

Find out more  
[www.philips.com/innospirego](http://www.philips.com/innospirego)

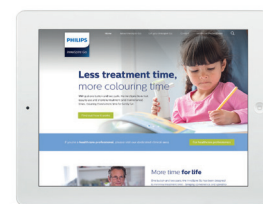

← Actual size →

\* Versus predecessor Aeroneb Go.

\*\* Using 2.5ml salbutamol.

\*\*\* A mask is included for children that find it difficult to use the mouthpiece.