

Consumable parts

Please refer to *Section 4, "Cleaning and maintenance"*, of your myAIRVO 2 User Manual for information on consumable parts. myAIRVO 2 has five parts which you will need to replace regularly to keep your device as clean as possible. The times recommended below assume your myAIRVO 2 is only being used by one person, and rely on the correct cleaning procedures being followed.

30 Optiflow™+ tracheostomy interface

Replace every 30 days.
Product code: OPT970E

Reorder code:

60 Heated breathing tube

Replace every 60 days.
Product code: 900PT560E (top right)
900PT500E (bottom right)

Reorder code:

60 MR290 auto-fill water chamber

Replace every 60 days.
Product code: 900PT290E

Reorder code:

60 Water bag

Replace every 60 days.
Product code: 900PT401 (2-pack)

Reorder code:

90 Air filter

Replace every 90 days or 1000 hours.
Product code: 900PT913 (2-pack)

Reorder code:

Some products may not be available in your country. Please contact your local Fisher and Paykel Healthcare representative.

Contact details for support and ordering

www.fphcare.com/myairvo

Important note: This is a Quick Start Guide only and should not be considered in any way a replacement for your myAIRVO 2 User Manual. Before using your myAIRVO 2 device you must read your myAIRVO 2 User Manual to understand the device's intended use, warnings and cautions. myAIRVO 2 is not intended for life support.

Setup steps for your myAIRVO 2 with the MR290 auto-fill water chamber

Please refer to *Section 2, "Setting up myAIRVO 2"*, and *Section 3, "Using myAIRVO 2"*, of your myAIRVO 2 User Manual for information on setup.

<p>1</p> <p>Place your myAIRVO 2 on a low shelf or near the floor beside your bed. Make sure it is placed below your head height and flat.</p>	<p>2</p> <p>a. Remove the blue port tabs and water supply tube holder. b. Install the provided elbow adapter to the MR290 water chamber. c. Clip the water supply tube into position.</p>
<p>3</p> <p>Fit the water chamber onto your device by depressing the finger guard, pushing forward until the finger guard clicks into place.</p>	<p>4</p> <p>a. Hang water bag from the home stand. b. Push water supply tube spike into the fitting at the bottom of the water bag. Open vent cap. c. Water will enter MR290 water chamber and be contained below the max water level line.</p>
<p>5</p> <p>Attach the heated breathing tube to your myAIRVO 2.</p>	<p>6</p> <p>Plug myAIRVO 2's power cord into the mains/utility power socket and into the rear of the unit.</p>
<p>7</p> <p>Switch on myAIRVO 2 by pressing and holding the On/Off button for 5 seconds.</p>	<p>8</p> <p>myAIRVO 2 will begin to warm up, displaying a warm-up symbol on the screen.</p>
<p>9</p> <p>When myAIRVO 2 is ready for you to use, a tick symbol will be displayed.</p>	<p>10</p> <p>Connect your interface and begin therapy. During use, if excess condensate forms in the tube, refer to <i>Condensate Management</i> in <i>Section 2, "Setting up myAIRVO 2"</i>.</p>

REF 614129 REV A November 2017 © Fisher & Paykel Healthcare Limited



Understanding myAIRVO 2 buttons

Please refer to Section 3, "Using myAIRVO 2", of your myAIRVO 2 User Manual for information on how to operate your myAIRVO 2.

	<p>On/Off button The On/Off button turns myAIRVO 2 on when held for 5 seconds. After use, press for 3 seconds to begin Drying Mode. Your device will turn off automatically after 99 minutes.</p>		<p>Audio Pause button The Audio Pause button allows you to temporarily silence audio alarms.</p>
	<p>Mode button The Mode button allows you to cycle between modes:</p> <ul style="list-style-type: none"> • Temperature • Flow Rate • Day/Night • Compliance 		<p>Up and Down buttons The Up and Down buttons allow you to change settings within the following modes:</p> <ul style="list-style-type: none"> • Temperature • Flow Rate • Day/Night


Understanding myAIRVO 2 modes and settings in default mode

Please refer to Section 3, "Using myAIRVO 2", of your myAIRVO 2 User Manual for information on how to operate your myAIRVO 2. Your myAIRVO 2 will remember your settings when switched off.

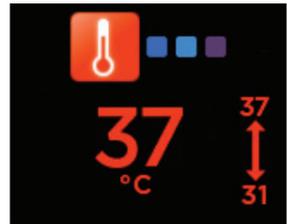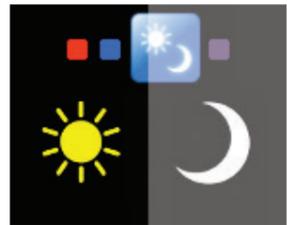	<p>Temperature Mode Please refer to Target Dew-Point Temperature in Section 3, "Using myAIRVO 2".</p>		<p>Flow Rate Mode Use the Up and Down buttons to select your flow rate from 10-60 L/min. Please see your myAIRVO 2 User Manual for information about flow ranges for your interface.</p>
	<p>Day/Night Mode Use the Up and Down buttons to cycle between Day Mode and Night Mode. In Night Mode some audio signals will become quieter and your display screen will dim.</p>		<p>Compliance Mode Compliance Mode automatically cycles between device usage statistics:</p> <ul style="list-style-type: none"> • Total hours used • Hours per day • Checksum (for your clinician)

myAIRVO 2 visual and auditory alarms

Please refer to Section 3, "Using myAIRVO 2", of your myAIRVO 2 User Manual for information on alarms. Your myAIRVO 2 has alarms to warn you about interruptions to your treatment. You can use the Audio Pause button to silence the alarm for 115 seconds.
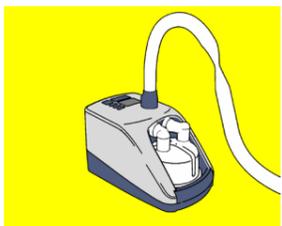
<p>Audio Pause button</p>	 <p>Check for blockages</p>	
----------------------------------	---	---

Cleaning myAIRVO 2


Please refer to Section 3, "Using myAIRVO 2", and Section 4, "Cleaning and maintenance", of your myAIRVO 2 User Manual for information on how to clean myAIRVO 2.

Daily cleaning

1



2

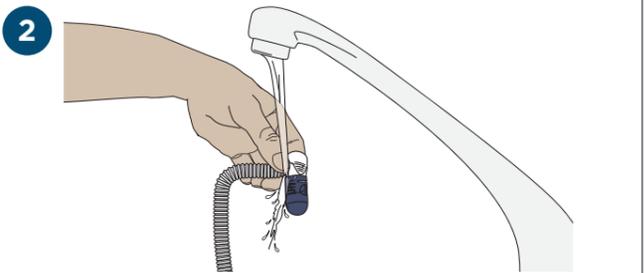


After use, remove your interface. To begin the 99-minute Drying Mode, press and hold the On/Off button for 3 seconds. To avoid burns, do not wear the interface during Drying Mode and do not remove the water chamber until Drying Mode is complete.


Rinse the interface in drinking-quality water then reconnect to the heated breathing tube whilst still in Drying Mode to dry the interface.

Weekly cleaning


1

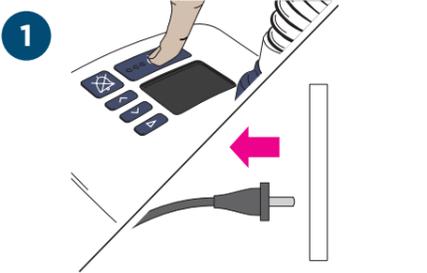


2

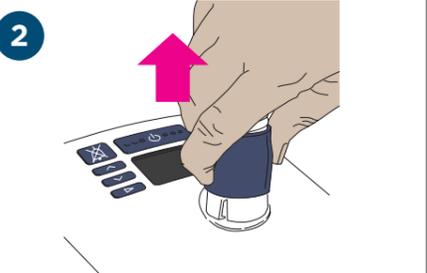

3

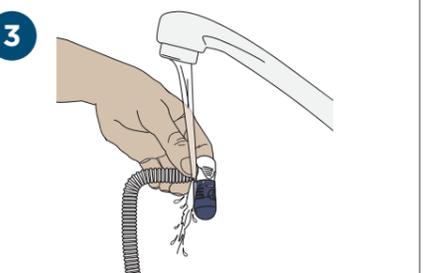

4

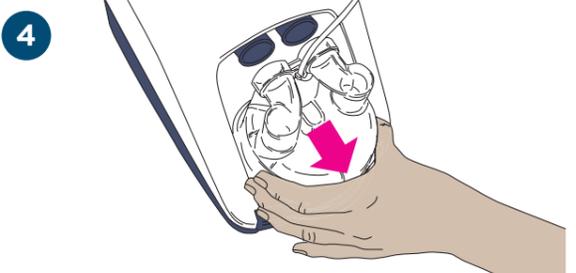

5



When cleaning is complete, set up your myAIRVO 2 using the setup steps previously described.