



**ResMed**


 **Respshop**  
better sleep, better life.  
**Click Here To Buy**

**AirTouch™ F20**




## Our softest mask ever

When it comes to effective sleep therapy, comfort is key. That's why the AirTouch F20 features the ResMed UltraSoft™ memory foam cushion, a breakthrough in CPAP mask comfort technology. This innovative new cushion conforms to the unique contours of each face for a truly personalized fit. It is also designed to deliver a light, breathable seal across a wide range of therapy pressures, fitting 98% of patients.<sup>1</sup>

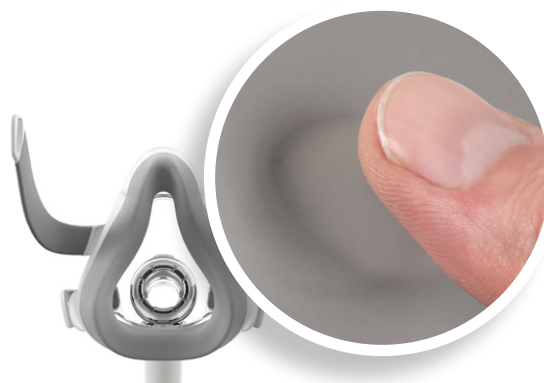

[ResMed.com/AirTouch20](https://ResMed.com/AirTouch20)

Supple memory foam adapts to the curves and contours of each face, creating a soft, personalized fit designed to reduce facial markings and skin irritation.

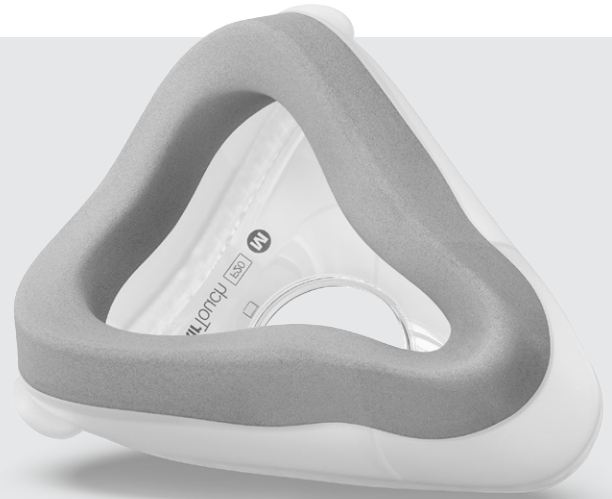

## AirTouch F20 highlights

### Plush headgear

specifically designed for extra softness and stability.

### UltraSoft

**memory foam cushion** is the softest cushion we've ever created, designed for exceptional comfort.

### Magnetic clips

quickly and easily guide the headgear to the frame.

### Soft and flexible frame

sits underneath the eyes and adapts to different facial types, while integrated padding adds softness. Compatible with InfinitySeal™ cushions.

### Diffused venting

evenly disperses air and dampens air flow.

### Quick-release elbow

makes it simple to disconnect from tubing when getting up during the night.





**AirTouch F20 for Her**  
Designed specifically for women.





## AirTouch F20 components

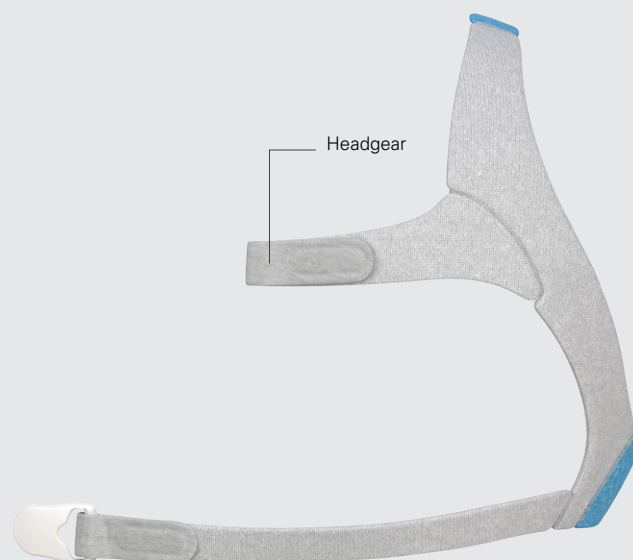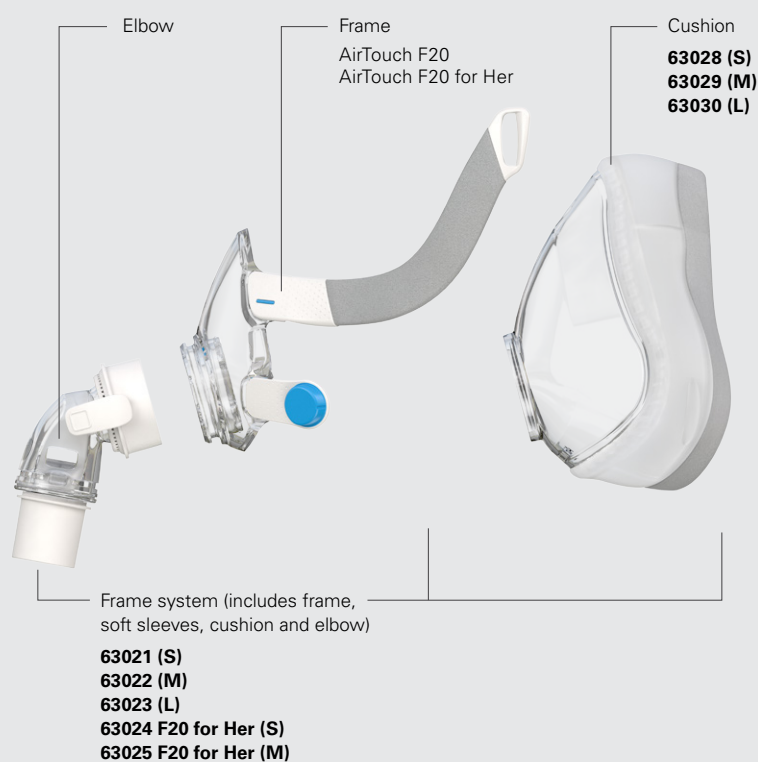

AirTouch F20 headgear  
with headgear clips  
**63470 (S)**  
**63471 (Std)**  
**63472 (L)**



AirTouch F20 for Her headgear  
with headgear clips  
**63473 (S)**

## Product codes

### Complete mask systems *(Includes frame system + headgear)*

AirTouch F20: Small	<b>63000</b>
AirTouch F20: Medium	<b>63001</b>
AirTouch F20: Large	<b>63002</b>
AirTouch F20 for Her: Small	<b>63003</b>
AirTouch F20 for Her: Medium	<b>63004</b>

**AirTouch™** F20



**AirTouch™ F20**



## Fitting guide for AirTouch F20



- Start by ensuring the magnets on the lower straps are unclipped.
- If they are clipped, detach them by twisting and pulling them away from the frame.



- Pull the headgear over the patient's head, making sure the logo is facing outwards.



- Position the lower straps so they both sit below the patient's ears.
- Bring each of the magnets up to meet their corresponding clips on the frame.



- Adjust the fastening tabs on the upper straps until the mask feels even and stable.
- Reattach the fastening tabs.



- Repeat step 4 with the lower straps.
- Reattach the fastening tabs.



- Connect the device's air tube to the elbow, then click the elbow into the front of the mask.
- Start the device.



### Additional options

For patients with longer hair, we recommend feeding their hair through the back of the headgear. If tied in a pony tail, it should sit high enough to fit comfortably through the back of the headgear.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

1 ResMed AirTouch F20 internal USA fitting study of new and existing CPAP patients, conducted April 2016 and August 2016.

Distributed by ResMed Corp, 9001 Spectrum Center Boulevard, San Diego, CA 92123 USA +1 858 836 5000 or 1 800 424 0737 (toll free). See ResMed.com for other ResMed locations worldwide. AirTouch, InfinitySeal and UltraSoft are trademarks and/or registered trademarks of the ResMed family of companies. Specifications may change without notice. © 2017 ResMed Ltd. 1019461/3 2017-01