RespShop

Resmed

AirTouch™ F40

AirFit™ F40

℞ Only


Clinical use only

Disinfection guide

This guide isintended for multi-patientre-useof only the SleepLab Mask (SLM)variants ofthe AirTouch F40 and AirFit F40 masks in a sleep lab, clinic or hospital. If you use themaskasa singleuser inthe home,referto the UserGuideforcleaning instructions.The processor isresponsiblefor ensuring that reprocessing is completed in accordance with Resmed's validated procedures. Potable water with hardness ≤150 mg/Lisrecommended for all cleaning, rinsing and disinfection steps.

	Cleaning Alconox™ (Mild alkaline, anionic detergent) Thermal Disinfection Manual 167°F (75°C) for 30 minutes 194-199°F (90-93°C) for 1 to 10 minutes	
Mask component[1]		Validated number of cycles
AirFit F40 cushion	✓	30
AirTouch F40 cushion	✓	
Frame	✓	
Multi-hole vent ring	✓	
Elbow and short tube	✓	
Headgear	✓	
QuietAir™ vent ring	-	Not intended for multi-patient use

[1] May not be available in all regions. For full details regarding the correct use of these masks, please refer to the User Guide. For a list of available replacement parts, check the Product Guide on Resmed.com

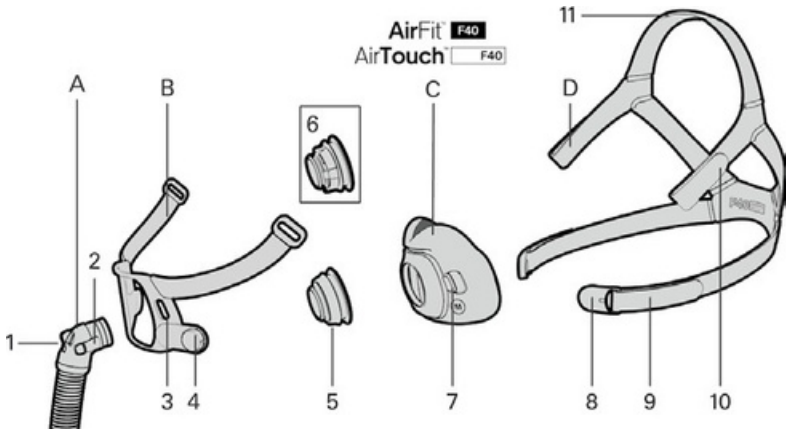



A	Elbow and short tube	5	QuietAir vent ring
B	Frame	6	Multi-hole vent
C	Cushion	7	ring Cushion stud
D	Headgear	8	Magnetic clip
1	Anti-asphyxia valve (AAV)	9	Lower headgear
2	Side button	1	strap Upper
3	Frame stud hole	0	headgear strap
4	Frame magnet	1	Top headgear
		1	strap

Manual cleaning and disinfection procedure

1. Disassembly	Disassemble themask accordingtotheinstructionsinthe User Guide.
2. Cleaning	1. Make a solution of a mild alkaline, anionic detergent and water as per manufacturer's instructions. Resmed has validated: • **Alconox** 1% concentration (10 g/L) 86-131°F (30-55°C) 2. Soak all mask components for 5 to 10 minutes. Ensure there are no air bubbles. 3. After the soaking time, while keeping the components in the solution, clean the inside and outside with a soft bristle brush as follows. A soft bristle tube/bottle brush is require the multi-hole elbow and short tube. Pay particular attention to the crevices and cavities. • **AirFit F40 cushion / AirTouch F40 cushion / frame / headgear** 1 minute brushing • **Multi-hole vent ring** 1 minute brushing, followed by 1 minute agitating in the solution. Pay particular attention to the vent holes. • **Elbow and short tube** 30 seconds brushing including minimum 15s of brushing inside of the short tube 4. For the AirTouch F40, rinse the cushion by shaking vigorously in 5 liters of water at 68-140°F (20-60°C) for 1 minute. 5. Rinse the components as follows under running water at 68-140°F (20-60°C) directly contacting all surfaces. • **AirFit F40 cushion / multi-hole vent ring / elbow and short tube** 30-60 seconds rinsing • **AirTouch F40 cushion** Repeat rinsing for 1 minute, while squeezing fabric part repeatedly • **Frame** 30-60 seconds rinsing, while squeezing the fabric component repeatedly and thoroughly • **Headgear** 90 seconds rinsing, while squeezing repeatedly and thoroughly
3. Inspection and drying	1. Inspect and if required, repeat all steps in section [2. Cleaning] until visually clean. 2. Squeeze the fabric components to remove excess water. 3. Allow the mask components to air dry out of direct sunlight. Drying is not required if thermal disinfection is continued immediately.
4. Disinfection	1. Fully immerse the mask components in a hot water bath using one of the following temperature-time combinations. Agitate the components in the water bath to ensure no air • 167°F (75°C) for 30 minutes, OR • 194-199°F (90-93°C) for 1 to 10 minutes 2. Squeeze the fabric components to remove excess water. 3. Allow the components to air dry out of direct sunlight.
5. Inspection	Perform a visual inspection of each mask component. If any visible deterioration of a mask component is apparent (cracking, crazing, tears, pilling etc.), the mask component shou
6. Reassembly	replaced. Slight discolouration of the silicone and fabric components may occur and is acceptable. Disinfecting fabric components together may result in dye transfer stains and ar
7. Storage and transportation	Reassemble the mask according to the instructions in the User Guide. For storage and transportation details, check the User Guide.

⚠ *GENERAL WARNINGS AND CAUTIONS*

- *Alwaysfollowcleaning and disinfection instructions. Some cleaning products may damage the mask, its parts and their function, or leave harmful residual vapors. Any deviations from the procedures listed in this guide, including exceeding the specified number of maximum cycles, can have an adverse effect on the mask and consequently the safety or the quality of therapy.*
- *Failure to thoroughly clean and disinfect the mask component as indicated may result in inadequate disinfection.*
- *When using detergents, disinfectants or equipment, always follow the manufacturer's instructions. In the event of conflict, this guide takes precedence.*
- *Always follow safe operating practices, including the use of appropriate personal protective equipment (PPE), as required. Refer to the manufacturer's safety instructions for more details.*

MANUFACTURER
ResMed Pty Ltd
1 Elizabeth Macarthur Drive
Bella Vista NSW 2153 Australia

See Resmed.com for other Resmed locations worldwide. AirFit, AirTouch and QuietAir are trademarks and/or registered trademarks of the ResMed family of companies. Alconox is a trademark of A other intellectual property information, see Resmed.com/ip. © 2026 ResMed. 608740/2 2026-05