




Performance  
Made Compact

# The next generation of Dynamic Support Technology™

The F&P Evora™ Full Mask is uniquely designed to optimize performance throughout the night.

Dynamic Support Technology combines Stability Wings with the Floating Seal to allow for freedom of movement while keeping the mask in place.

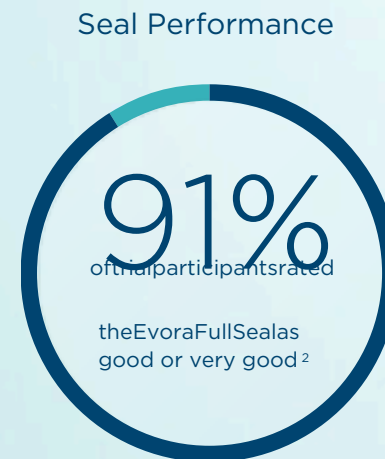


## Secure Comfort

The Floating Seal accommodates for movement while the soft silicone gently wraps under the nose for a comfortable night's sleep.




Comfort

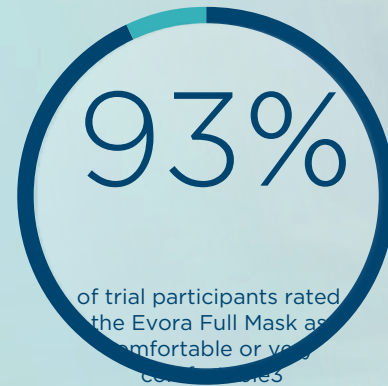




Removable Frame Clip  
Enables the headgear settings to stay in place during fitting and cleaning

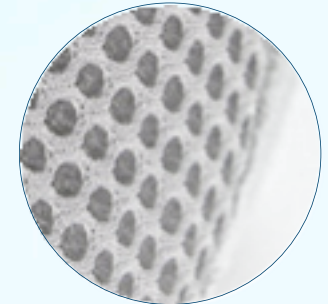

Headgear with VentiCool™ Technology  
Breathable headgear designed to allow for heat and moisture to escape

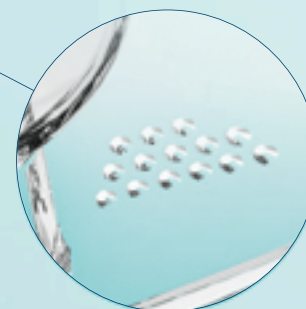

Exhaust Holes  
Designed for minimal noise



Sizing Options



Extra Small



Small to Medium



Large



Standard



Extra Large

Fitting Summary



Guide the headgear over your head



Hook both the headgear clips onto the frame



Adjust top, then bottom straps

For further instructions, please refer to Use and Care Guide.

Product Codes


Code	Description
EVF1XA	Evora Full Mask Extra Small A Model
EVF1MA	Evora Full Mask Small-Medium A Model
EVF1LA	Evora Full Mask Large A Model Evora
EVF1XMLA	Full Mask Fitpack A Model

Spare Parts

Universal frame fits all seal sizes	
Code	Description
400EVF114	Evora Full Seal - Extra Small Spare
400EVF115	Evora Full Seal - Small to Medium Spare
400EVF116	Evora Full Seal - Large Spare
400EVF121	Evora Full Headgear - Standard Spare
400EVF122	Evora Full Headgear - Extra Large Spare
400EVF131	Evora Full Mask - No Headgear, Extra Small Seal
400EVF132	Evora Full Mask - No Headgear, Small to Medium Seal
400EVF133	Evora Full Mask - No Headgear, Large Seal
400EVF141	Evora Full Headgear Clips & Frame Clip Spare
400EVF151	Evora Full Tube Spare
400EVF161	Evora Full Frame Spare
400FPH341	F&P Swivel White Spare

F&P Evora Full is available on the F&P myMask™ App

Instructs patients on mask use through step-by-step fitting, fine-tuning and cleaning videos. Download the F&P myMask App or scan the QR code to download directly.

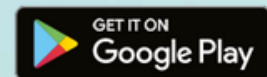

1. 48 OUT OF 44 PARTICIPANTS RATED THE EVORA FULL STABILITY AS EITHER VERY STABLE OR STABLE. INTERNAL VALIDATION STUDY CONDUCTED WITH 44 PARTICIPANTS IN NEW ZEALAND 2020. EVORA FULL MASK AS VERY GOOD OR GOOD. INTERNAL VALIDATION STUDY CONDUCTED WITH 44 PARTICIPANTS IN NEW ZEALAND 2020. 3. 41 OUT OF 44 RATED THE EVORA FULL MASK AS VERY GOOD OR GOOD. INTERNAL VALIDATION STUDY CONDUCTED WITH 44 PARTICIPANTS IN NEW ZEALAND 2020.

APPLE, THE APPLE LOGO AND IPHONE ARE TRADEMARKS OF APPLE INC., REGISTERED IN THE U.S. AND OTHER COUNTRIES. APP STORE IS A SERVICE MARK OF APPLE INC., REGISTERED IN THE U.S. AND OTHER COUNTRIES. GOOGLE PLAY AND GOOGLE PLAY LOGO ARE TRADEMARKS OF GOOGLE LLC. F&P, EVORA, MYMASK, DYNAMIC SUPPORT TECHNOLOGY AND VENTICOOL ARE TRADEMARKS OF FISHER & PAYKEL HEALTHCARE LIMITED. FOR PATENT INFORMATION, PLEASE SEE WWW.FPHCARE.COM/IP. REF 625824 REV A 2022-03 © 2022 FISHER & PAYKEL HEALTHCARE LIMITED