


SPIemrfpolrym Fsi tAs u+t oPmerafotircmalsy - in Tune with You

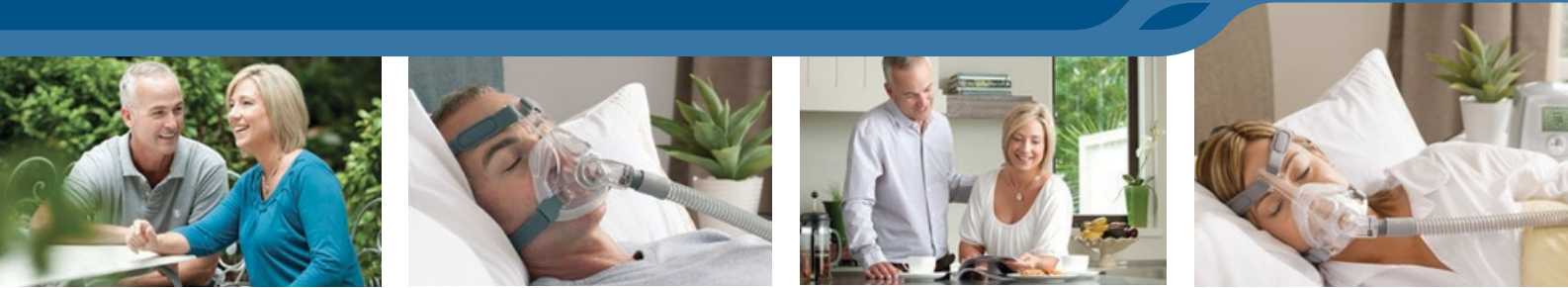

**F&P SIMPLUS FULL FACE MASK**

*When engineering the F&P Simplus our designers set out to create a mask that revolutionized Full Face comfort, seal and easy use.*

*The F&P Simplus is the result of intensive research and development, many prototypes and thousands of design hours. We consulted patients, equipment providers and clinicians from around the world to develop a mask that simply fits + performs.*

## F&P Simplus Features

The F&P Simplus incorporates three key components, the **RollFit™ Seal**, **ErgoForm™ Headgear** and **Easy Frame**, all designed to work in harmony. In combination, these components offer the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.



### 1. RollFit Seal

*As the name suggests, the one-piece seal “rolls” back and forth on the bridge of the nose to adjust automatically. In doing so, the RollFit Technology reduces pressure on the bridge of the nose.*

*The RollFit Seal incorporates Integrated Stability Panels to enhance the sealing performance and the Advanced Air Diffuser is designed for minimal noise and draft.*

### 2. ErgoForm Headgear

*This breathable, ergonomically designed headgear self-locates high on the rear of the head, allowing for unrestricted head movement (sideways, up and down).*

*Incorporating Stretch and Non-stretch Panels, the ErgoForm Headgear is designed to hold its shape and allow stretch where it is required, for optimal usability and performance.*

*Embedded VELCRO®\* Tabs with Easy-Clip Hooks enable fast fitting and removal.*

*The Adjustable Crown Strap enhances mask stability.*

### 3. Easy Frame

*This low-profile frame is stable, durable and small, and ensures a clear line of sight.*

*The one frame fits all three seal sizes and has an Easy-Clip Frame Attachment to aid assembly after cleaning.*

*The Ball-and-Socket Elbow rotates for freedom of movement while reducing drag from the CPAP tube.*

Mask Reorder Codes		Spare Parts Reorder Codes			
40047	F&P Simplus Mask – Small	400HC57	F&P Simplus Seal – Small	400HC58	F&P Simplus Headgear – Small
5	F&P Simplus Mask – Medium	9	F&P Simplus Seal – Medium	2	F&P Simplus Headgear – Medium/Large
40047	F&P Simplus Mask – Large	400HC58	F&P Simplus Seal – Large	400HC58	F&P Simplus Mask No H/G – Small**
6		0	F&P Simplus Swivel	3	F&P Simplus Mask No H/G – Medium**
40047		400HC58	F&P Simplus Headgear – Clips & Buckle	400HC58	F&P Simplus Mask No H/G – Large**
7				4	
		400HC57		400HC58	
		0		5	
		400HC58		400HC58	

\* VELCRO® is a registered trademark of Velcro Industries B.V.

\*\* Not available in all countries