

Now, for the first time, active, working, and traveling oxygen patients can truly enjoy the lifestyle they desire, with an unrestricted oxygen supply.

## Mobility Redefined

With AirSep's LifeStyle™ Portable Oxygen Concentrator, those who need the benefits of oxygen therapy can travel and vacation with ease, now including travel by air\*, or simply enjoy unplanned, overnight stays. Supplemental oxygen users can go wherever adventure calls – without the constant fear of running out of their liquid or cylinder supply. And it's now so easy to pack for a trip to just about anywhere – with one compact, lightweight (less than 10 lb/4.5 kg) unit that makes its own supply of oxygen, whether at home or away.

## No Refills, No Deliveries Ever Again!

The LifeStyle concentrator encapsulates the technology and dependability of a traditional concentrator with the economics of a pulse flow conserving device, providing your patients with five pulse settings, according to prescription. And what's more, you know the peace of mind free from the worry of deliveries and refills, particularly with a highly ambulatory or remotely-located patient.

With LifeStyle, your patient always has an endless supply of oxygen on hand because the device is truly a mini oxygen generating system, producing its own contents everywhere your patient goes.

## Multiple Power Choices

AirSep offers the most flexible oxygen supply in today's market, powered by normal household AC electric, DC automobile power, or a rechargeable

battery that easily snaps in and out of the unit for quick change-outs.

In addition, an optional airplane DC adapter is available from AirSep for use onboard commercial aircraft that enable the use of portable oxygen concentrators through the plane's laptop connection, as equipped.

## Marketing Advantages

Equipment providers and patients who use the portable oxygen concentrator on a daily basis easily



see two primary benefits and a winning combination in LifeStyle: economics and convenience. You can fully promote mobility while eliminating the costly logistic issues of supplying compressed gas or liquid portable oxygen systems. And patients appreciate all of their oxygen needs met at a fraction of the size and power consumption of traditional concentrators.

The priceless benefits of increased ambulation, essential to a patient's health and well-being, is now within reach.

Endless lifestyle possibilities for the patient and the provider...

*\*As of August 2005, the US Federal Aviation Administration (FAA) approved LifeStyle for in-flight use (FAA2004-18596, SFAR No. 106), according to a physician's prescription, thereby enabling individual US airlines to determine and establish portable oxygen concentrator policy.*

*In addition, airlines based in other countries also may enable oxygen patients to fly with LifeStyle, depending on the air carriers' own determination on how they are best equipped to offer this convenience to traveling oxygen patients.*

*As always, check directly with the medical oxygen desk of any airlines your patients may wish to travel with prior to their trip to ensure that the policy of the particular airline allows both carry on and in-flight use of LifeStyle Portable Oxygen Concentrators.*



## For Active and Traveling Patients...

LifeStyle's multiple power choices enable total flexibility within one versatile, easy-to-use, easy-to-handle unit:

- ▶ AC adapter with universal electrical power supply (100-250 volt, 50/60 Hz)
- ▶ DC adapter (12-volt) for automobiles/motorized boats
- ▶ Optional DC adapter for onboard airline use
- ▶ Rechargeable, internal battery

Simply change the detachable power cord to LifeStyle's universal AC power supply, as necessary, and you have a unit that can travel with your patients — around the corner or around the globe.

## ...For Working Patients

A single, portable, contemporary-styled unit that makes and supplies all of a patient's oxygen needs — all day long — is ideal for use in most full-time or structured work situations.


Even part-time or semi-retired individuals on oxygen therapy can remain gainfully employed or choose to do volunteer activities with a little help from their LifeStyle. LifeStyle is suitable for working at open-air markets, sitting behind the desk in a busy office, cashiering at a store's checkout counter, and a host of other work activities — both indoors and out.



## ...Or for Leisure Pursuits

Whether going to the beauty parlor, a neighbor's home, a sporting or entertainment event, relaxing with a favorite hobby, or just enjoying the convenience of portable

*LifeStyle*™  
by AirSep®

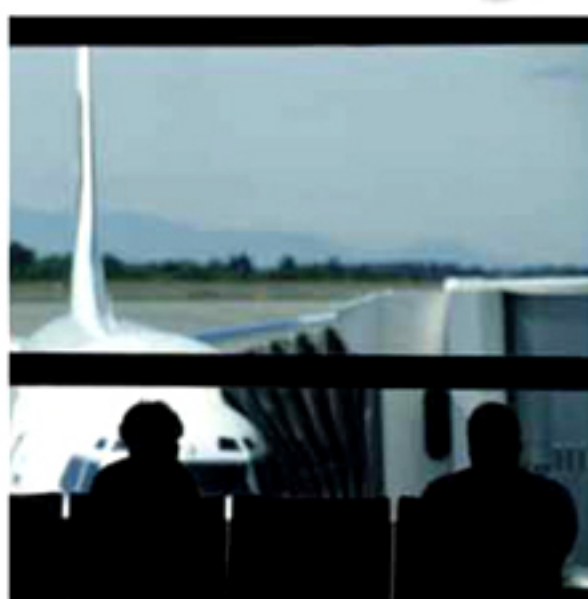

oxygen without the fuss of managing contents and rationing supplies, LifeStyle is the answer for today's savvy oxygen user.

The long-awaited flight across the country — or even across the world — to see distant relatives or dear friends — may be closer at hand with the help of LifeStyle.

## Convenience at Home...

Even moving about the house has never been easier — with the oxygen concentrator that patients can actually carry over their shoulder while walking. Climbing up and downstairs while LifeStyle makes it own oxygen — is total confidence for the O<sub>2</sub> patient!





## SPECIFICATIONS

<b>OXYGEN CONCENTRATION*:</b>	1-5 pulse settings; equivalent to a continuous flow of 90% oxygen at $\pm 3$
<b>DIMENSIONS:</b>	5.5 in. H x 7.25 in. W x 16.31 in. L (13.9 cm H x 18.4 cm W x 41.4 cm L)
<b>WEIGHT:</b>	9.75 lb (4.4 kg)
<b>POWER:</b>	LifeStyle unit – 13 VDC, 42 Watts Power supply – 100-250 VAC, 50/60 Hz, 2.0 amp maximum Battery – 12 VDC, nominal
<b>BATTERY DURATION:</b>	50 minutes nominal for fully-charged battery; up to 3-1/4 hours with 3 batteries in optional PowerPack
<b>BATTERY CHARGE TIME:</b>	2 to 2-1/2 hours; PowerPack with up to 3 batteries – 4 hours
<b>AUDIBLE ALARMS AND VISUAL INDICATORS:</b>	Start-up – Audible/visual (green) Pulse flow – Visual (green) Battery condition – Audible/visual (yellow) Cannula disconnect – Audible/visual (red) System overdraw – Audible/visual (red) High and low pressure – Audible/visual (red)
<b>WARRANTY:</b>	2-year parts on LifeStyle unit; 90 days on rechargeable battery; 2-year parts on PowerPack

\*Based on atmospheric pressure of 14.7 psia (101 kPa) at 70° F (21° C)

## INTERNATIONAL APPROVALS:

LifeStyle Unit –




- CE 0459
- CSA-certified
- FAA-approved
- RTCA DO160-D, Section 21, Category M-compliant
- AirSep Corporation is ISO 9001:2000-certified



**AirSep Corporation**  
**Medical Division**

401 Creekside Drive • Buffalo, New York 14228-2085 USA  
USA/Canada Toll-free: 800.874.0202  
Tel: 716.691.0202 • Fax: 716.691.4141  
E-mail: mpd@airsep.com • Internet: www.airsep.com

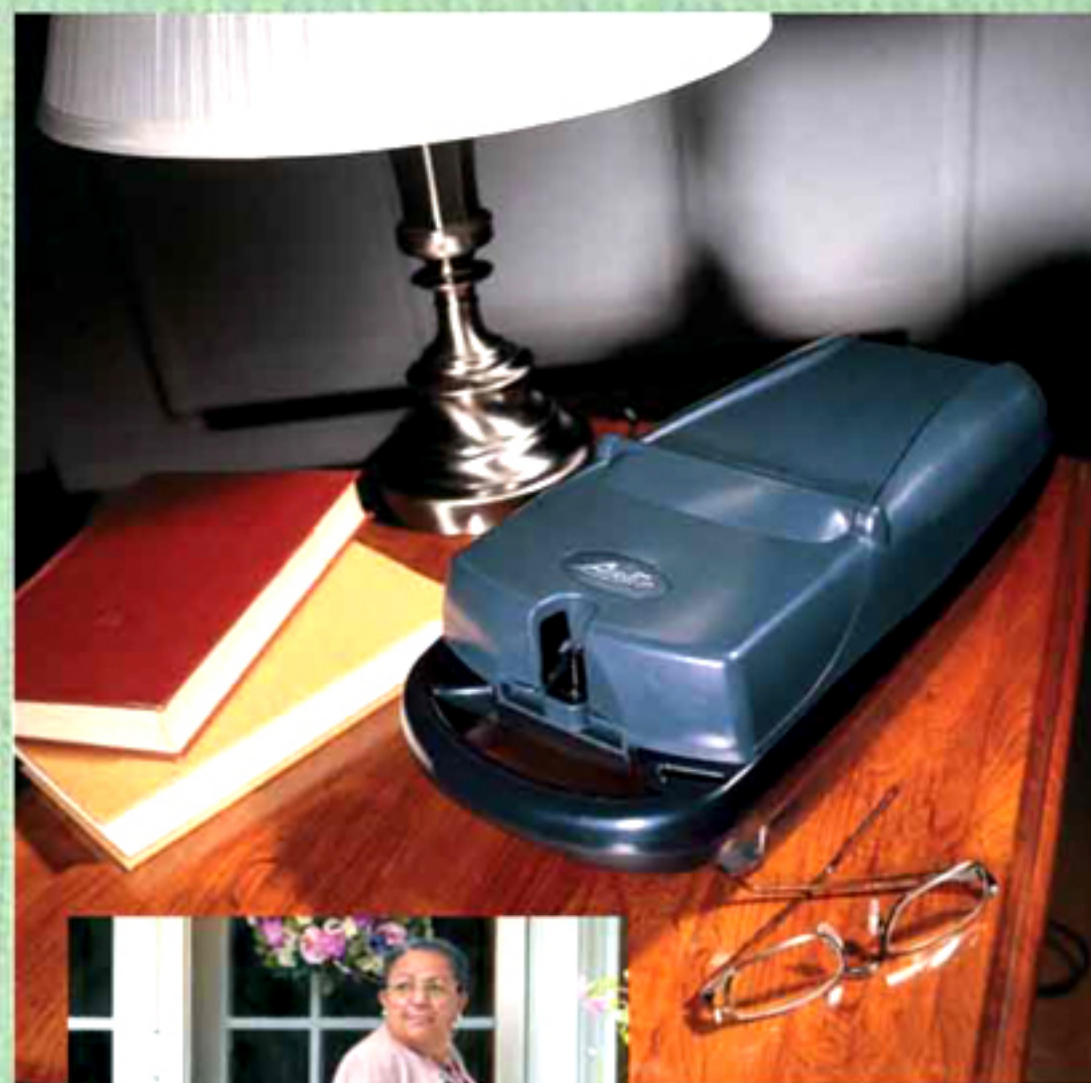

This optional lightweight, retractable cart makes transport of the LifeStyle concentrator and PowerPack (as shown) easy during extended use.

The optional PowerPack is used to recharge LifeStyle batteries and greatly increases ambulation time.

## ORDERING INFORMATION

### PART NO. DESCRIPTION

**AS081-1 LifeStyle Portable Oxygen Concentrator:**  
Includes unit, comfortable shoulder strap/carrying bag, rechargeable battery, AC power supply (PW005-1) with detachable power cord, DC automobile adapter (EL254-1), and patient manual.

### Additional travel accessories:

<b>BT007-1</b>	<b>Rechargeable battery</b>
<b>BT011-1</b>	<b>PowerPack battery charger (120 Volt)</b>
<b>BT011-3</b>	<b>PowerPack battery charger (220-240 Volt)</b>
<b>CD023-1</b>	<b>Power cord (120 Volt)</b>
<b>CD023-2</b>	<b>Power cord (120 Volt hospital-grade)</b>
<b>CD017-2</b>	<b>Power cord (220-240 Volt)</b>
<b>CD025-1</b>	<b>Power cord (220-240 Volt for Australia)</b>
<b>CD030-1</b>	<b>Airline-compatible DC adapter</b>
<b>MI210-1</b>	<b>Retractable cart</b>