



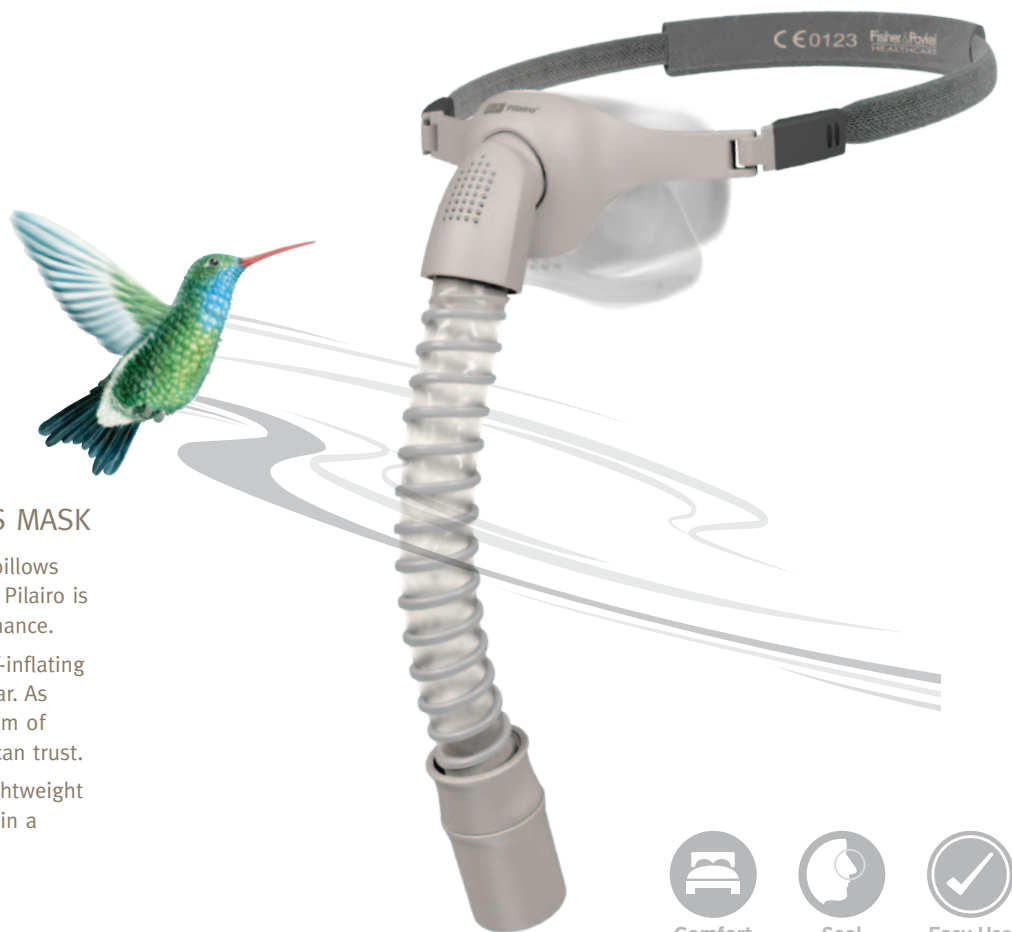

Sleep freely on pillows of air

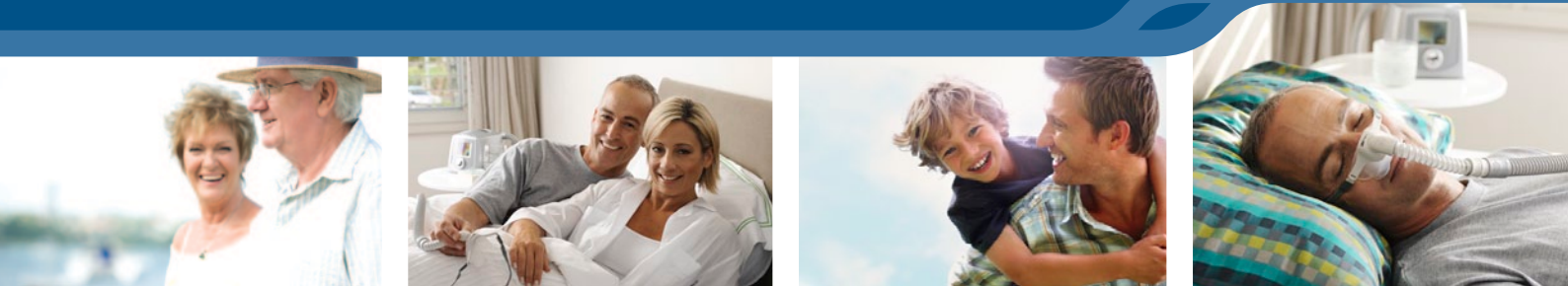

### F&P PILAIRO NASAL PILLOWS MASK

Weighing in as the lightest F&P nasal pillows mask (1.83 ounces/52 grams), the F&P Pilairo is light on the patient and big on performance.

The clever design integrates a new self-inflating AirPillow™ Seal and minimalist headgear. As a result, the patient experiences freedom of movement coupled with stability they can trust.

The revolutionary AirPillow Seal and lightweight comfort of the F&P Pilairo truly puts it in a league of its own.



Comfort



Seal



Easy Use

Every detail of the F&P Pilairo has been carefully researched, developed and engineered to improve patient SleepLife\*. The F&P Pilairo represents minimalist, highly effective nasal pillows design at its best and sets new benchmarks for comfort, seal and easy use.



### Feel the new AirPillow Seal.

- The micro-fine contoured silicone used in the F&P Pilairo's AirPillow Seal is (0.04 inches/0.12 mm) at its thinnest point and only (0.12 inches/3.0 mm) at its thickest.
- Soft, comfortable and highly effective, this fine seal is ergonomically designed to self-inflate with (Continuous Positive Airway Pressure) CPAP flow, to gently envelop the nose.



### Experience total comfort on pillows of air.

- The F&P Pilairo makes the need for complicated adjustments and complex headgear redundant.
- Its minimalist design is stable, low profile and unobtrusive.
- No more messy straps, adjustments or facial indents.

### Experience easier therapy.

- One convenient size.
- Only three simple parts to clean.
- Quick and easy to fit.



Comfort



Seal



Easy Use

### A little bird told us to make it light, small, fast and clever.

The design of the F&P Pilairo was inspired by the aerodynamic flight of the world's lightest little bird - the hummingbird. The hummingbird can rotate its fine wings in a circle; it is the only bird able to fly forwards, backwards, up, down and sideways and can hover in midair. This drove the design of a mask seal which self-inflates around the nares in a 'hover-like way' while providing total freedom of movement.





This is our lightest mask, the F&P Pilairo weighing in at 1.83 ounces/ 52 grams.

**1. AirPillow Seal**

- The AirPillow Seal is 3D contoured with fine Variable Thickness Silicone (VTST™).
- It uses CPAP pressure to self-inflate to create an effective yet gentle seal into and around the nose. The AirPillow Seal can move and flex in any direction without compromising the seal.

**2. Stabilizers**

- The stabilizers provide the support frame for the AirPillow Seal and assist with stability.

**3. Stretchwise™ Headgear**

- Made from a soft latex-free elasticized thread with intelligent-stretch properties, the headgear automatically customizes to the wearer.
- No manual adjustments are necessary.
- One convenient size.

**4. Flexible Tube**

- The delivery tube is lightweight, very flexible and designed to minimize pressure or pull on the face.
- The Ball and Socket Elbow provides greater flexibility for tube positioning.

**5. Quiet Air Diffuser**

- Efficiently designed to ensure the flow is directed away from the patient and their bed partner for a quiet sleep.

FITTING THE F&P PILAIRO IS EASY. FOLLOW THE SIMPLE SEQUENCE BELOW FOR BEST RESULTS:



**1. Spread the Wings**  
Use your fingers to pull apart the pillows.



**2. Position the Mask**  
Place the pillows on your nose until you feel a comfortable seal.



**3. Locate Headgear**  
Pull on the auto-adjusting, one-size stretchwise headgear.



**4. Test the Wings**  
Give the wings (pillows) a little wiggle for comfort.



**5. Turn on CPAP device**  
The F&P Pilairo's AirPillow Seal will self-inflate and seemingly hover in and around your nose.

For more information please contact  
your local Fisher & Paykel Healthcare representative

**Manufacturer** 

Fisher & Paykel Healthcare Ltd  
15 Maurice Paykel Place  
East Tamaki, Auckland 2013

PO Box 14 348 Panmure  
Auckland 1741  
New Zealand

Tel: +64 9 574 0100  
Fax: +64 9 574 0158  
Email: [info@fphcare.co.nz](mailto:info@fphcare.co.nz)  
Web: [www.fphcare.com](http://www.fphcare.com)

**Australia**

Fisher & Paykel Healthcare Pty Limited  
36-40 New Street, PO Box 167  
Ringwood, Melbourne  
Victoria 3134, Australia

Tel: +61 3 9879 5022  
Fax: +61 3 9879 5232

**Austria**

Tel: 0800 29 31 23  
Fax: 0800 29 31 22

**Benelux**

Tel: +31 40 216 3555  
Fax: +31 40 216 3554

**China**

Tel: +86 20 3205 3486  
Fax: +86 20 3205 2132

**France**

Tel: +33 1 6446 5201  
Fax: +33 1 6446 5221

**Germany**

Tel: +49 7181 98599 0  
Fax: +49 7181 98599 66

**India**

Tel: +91 80 4284 4000  
Fax: +91 80 4123 6044

**Irish Republic**

Tel: 1800 409 011

**Italy**

Tel: +39 06 7839 2939  
Fax: +39 06 7814 7709

**Japan**

Tel: +81 3 3661 7205  
Fax: +81 3 3661 7206

**Northern Ireland**

Tel: 0800 132 189

**Spain**

Tel: +34 902 013 346  
Fax: +34 902 013 379

**Sweden**

Tel: +46 8 564 76 680  
Fax: +46 8 36 63 10

**Switzerland**

Tel: 0800 83 47 63  
Fax: 0800 83 47 54

**Taiwan**

Tel: +886 2 8751 1739  
Fax: +886 2 8751 5625

**Turkey**

Tel: +90 312 354 34 12  
Fax: +90 312 354 31 01

**UK** 

Fisher & Paykel Healthcare Ltd  
Unit 16, Cordwallis Park  
Clivemont Road, Maidenhead  
Berkshire SL6 7BU, UK

Tel: +44 1628 626 136  
Fax: +44 1628 626 146

**USA/CANADA**

Tel: +1 800 446 3908  
or +1 949 453 4000  
Fax: +1 949 453 4001

フィッシャー & パイケルヘルスケア株式会社  
〒103-0025 東京都中央区日本橋茅場町3-8-10 リベラ茅場町5F  
FAX. 03-3661-7206 ご注文書専用 FAX. 0800-100-7207  
TEL. 03-3661-7205 (代表)

**Contents**

- Complete mask with headgear assembled and ready for use
- One convenient size
- User-friendly tear off Fitting Guide

**Reorder Code**

400420 F&P Pilairo