

Sleep freely on pillows of air



## Freedom... with added features

- ⊕ **Q Cover** for quieter use
- ⊕ **Choice of headgear** for versatility
- ⊕ Still one **convenient mask size** with three simple parts to clean


It's just what you asked for – a quieter nasal pillows mask with the option of Adjustable or Stretchwise™ Headgear. This enhanced offering ensures quieter use, less draft and fits a greater range of head sizes.

Weighing in as one of the lightest CPAP nasal pillows masks (1.94 ounces/55 grams),\* the F&P Pilairo Q is quieter for the patient and bigger on performance. The clever design integrates the self-inflating AirPillow™ Seal and a choice between two headgears to accommodate more of your patients. So, offer your patients freedom of movement coupled with stability they can trust.

\*F&P Pilairo Q with Stretchwise Headgear



The F&P Pilairo Q has been developed to improve patient satisfaction. Designed for comfort, seal and easy use, the F&P Pilairo Q represents a highly effective nasal pillows mask with just three key components and two minimalist headgear options.

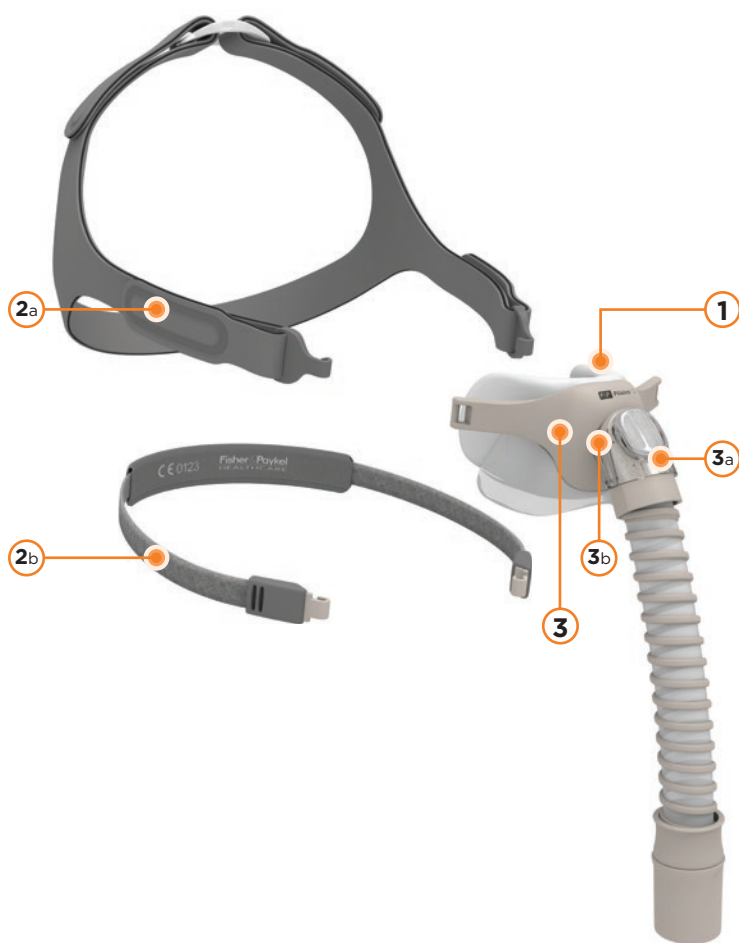

### Three Key Components:

#### 1. AirPillow Seal

The AirPillow Seal is 3D contoured with fine Variable Thickness Silicone (VTS™) and uses CPAP pressure to self-inflate to create an effective yet gentle double seal into and around the nose. The AirPillow Seal can move and flex in any direction without compromising the seal.

#### 2. Headgear Options: Adjustable or Stretchwise Headgear

##### a. Adjustable Headgear

This breathable, ergonomically designed headgear fits a greater range of head sizes and enables sizing adjustment to the patient's preference.

##### b. Stretchwise Headgear – an alternative option

Made from a soft, elasticized, latex-free thread with intelligent-stretch properties, this minimal headgear automatically adjusts to the wearer for a quick and easy fit.

#### 3. Easy Frame

This low-profile frame is stable, durable and small, and ensures a clear line of sight. The delivery tube is lightweight, very flexible and designed to minimize pressure or pull on the face. The stabilizers provide the support frame for the AirPillow Seal and assist with stability.

##### a. Q Cover

For quieter use<sup>1</sup> and minimal draft.

##### b. The Ball-and-Socket Elbow

Rotates for freedom of movement while reducing drag from the CPAP tube.

<sup>1</sup> The Sound Pressure Level of the mask with the Q Cover is 23.6 dBA, with uncertainty 2.5 dBA.



Reorder codes			
<b>400421</b>	F&P Pilairo Q – includes Adjustable Headgear and Stretchwise Headgear	<b>400HC565</b>	F&P Pilairo Q Nasal Pillows Mask (no headgear)
<b>400HC125</b>	F&P AirPillow Seal	<b>400HC229</b>	F&P Pilairo Q Elbow Cover & Diffuser (10 pack)
<b>400HC328</b>	F&P Pilairo Q Adjustable & Stretchwise Headgear (combo pack) – USA only	<b>400HC326</b>	F&P Pilairo Q Adjustable Headgear – Excludes USA
		<b>400HC588</b>	F&P Pilairo Q Headgear Clips & Buckle
		<b>400HC224</b>	Female Swivel
		<b>400HC322</b>	F&P Pilairo Q Stretchwise Headgear – Excludes USA