



Introducing the Next  
Generation BiLevel by 3B Medical.

## Luna G3 BPAP 25A

Integrated heated tubing and modem makes for a streamlined profile and a small footprint. The Luna G3 BPAP 25A is loaded with new features and settings to make bilevel therapy both more comfortable and convenient.



Luna G3 BPAP 25A

SKU: LG3700

GTIN: 00817783021159

 **3B™ Medical, Inc.**



**Respshop**  
better sleep, better life.

# Specifications

## DEVICE SIZE

Dimensions (L x W x H): 265 mm × 145 mm × 114 mm

Weight: 1.7 kg

Water capacity: To maximum fill line 360 mL

## PRODUCT USE, TRANSPORT AND STORAGE

Operation Transport and Storage

Temperature: 5°C to 35°C (41°F to 95°F) -25°C to 70°C (-13°F to 158°F)

Humidity: 15% to 93% Non-condensing 15% to 93% Non-condensing

Atmospheric Pressure: 760 ~ 1060 hPa

## MODE OF OPERATION

Continuous

## WORK MODE

CPAP, S, AutoS

## SD CARD

The SD card can record patient data and fault information

## AC POWER CONSUMPTION

100V - 240V ~, 50Hz / 60 Hz, 2 A

## DEVICE OFFER TO HEATED TUBING COMMUNICATIONS PORT

24V === 18W

## PRESSURE RANGE

4 to 25 cmH<sub>2</sub>O (in 0.5 cmH<sub>2</sub>O increments)

Under single fault conditions, ≤ 30 cmH<sub>2</sub>O for CPAP and AutoCPAP mode,

## RAMP

The ramp time ranges from 0 to 60 minutes.

## SOUND PRESSURE LEVEL

< 26 dB, when the device is working at the pressure of 10 cmH<sub>2</sub>O.

## THE FORM AND THE DIMENSIONS OF THE PATIENT CONNECTION PORT

The 22 mm conical air outlet complies with ISO 5356-1.

# Replacement Part Numbers

LG31020

Luna G3 Air Filter

LG34510

Luna G3 Water Chamber


LG3B

Luna G3 Carry Bag


LG3HT

Luna G3 Heated Tubing

# Ultra Modern Design with a User Friendly Interface

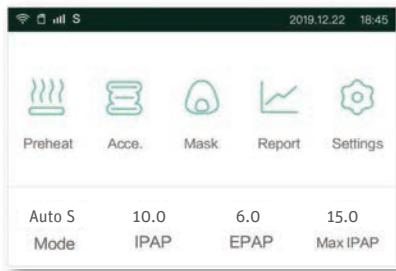

Easy access to frequently used features.

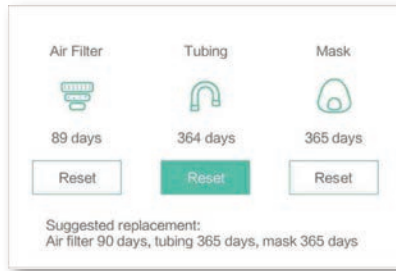



Accessory and supply reminders keep your treatment on track.



Clock, Temperature and Humidity Monitor on stand-by for convenience.

## Pre-Heat Feature

The G3 series has a pre-heat feature to warm your humidifier water up for ultimate comfort.

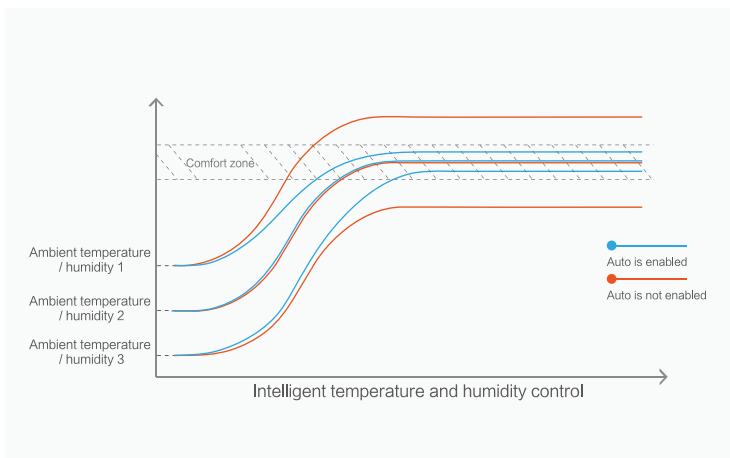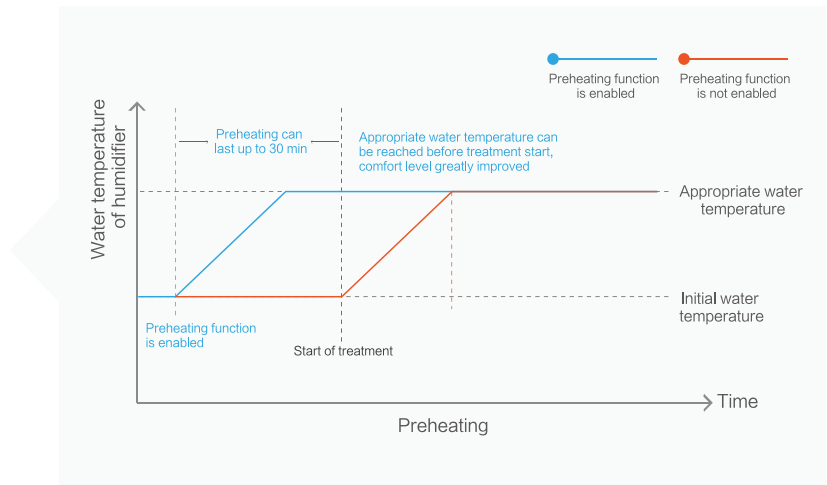


## Integrated Heated Tubing

intelligent Auto Humidity control allows G3 maximum moisture while ensuring rain out in the tubing is not an issue.

# Auto Ramp

New auto ramp feature senses when the user falls asleep and will adjust ramp time appropriately.

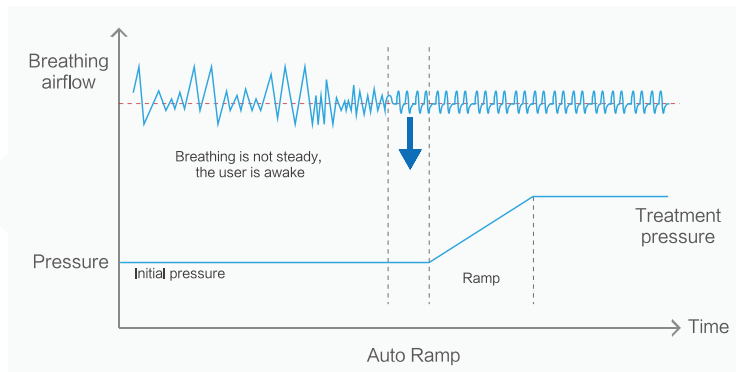

# Excellent Sound Performance

Multi-stage noise reduction technology is used to ensure sound performance is under 26 dB.

# Data Transfer

