



## A system you can rely on

### Philips Respironics InnoSpire Essence nebulizer compressor system

#### **A reliable choice**

When reliability is key, choose InnoSpire Essence. Essence, part of our InnoSpire family of compressor nebulizer systems, provides fast, efficient and reliable drug delivery to patients at home.

#### **Proven nebulizer technology**

Essence is matched with our SideStream disposable nebulizer to deliver a high quality aerosol in a short treatment time. With more than 20 years of aerosol therapy excellence, you can rest assured that patients will be provided with a high quality aerosol every time.


#### **Key features**

- Partnered with proven SideStream technology
- Short treatment times in 6 – 8 minutes\*
- Compact design



## Reliable aerosol delivery you can trust

Essence is combined with our highly efficient SideStream nebulizer. SideStream has been designed to draw in more air, via an active venturi system, resulting in faster drug delivery and shorter treatment times. It also has a unique diamond jet design that helps to minimize wear from frequent use and aids consistent delivery of a high quality aerosol.

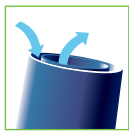

By utilizing a venturi, SideStream produces a flow rate of 16 l/min



Innovative jet design, producing consistently high quality aerosol

### Performance data

MMAD	2.9 µm
Fine particle fraction	77% below 5 µm
Treatment time (2.5 ml charge)	6 - 8 mins

### Product description (with voltage/frequency)

### Part number

Essence 230 V/50 Hz EU plug	1099967
Essence 220 V/60 Hz Korean plug	1099965
Essence 220 V/60 Hz UK plug	1099964
Essence 230 V/50 Hz UK plug	1112261
Essence 220 V/60 Hz US plug	1102174
SideStream reusable year pack	2010A
SideStream reusable kit with adult mask	1223A
SideStream reusable kit with pediatric mask	1224A
SideStream reusable kit with mouthpiece	1225A
Air inlet filters (pack of 4)	1102088
SideStream disposable one month pack (pack of 40)	2350
SideStream disposable with adult mask (pack of 50)	4446
SideStream disposable with pediatric mask (pack of 50)	4447
SideStream disposable with mouthpiece (pack of 50)	4448

### Technical specification

Electricity supply	230 V/50 Hz, 1 Amp	220 V/60 Hz, 1.2 Amp
Compressor max. pressure	317 kPa (46 psi)	296 kPa (43 psi)
Average flow rate	6 lpm @ 10 psig	7 lpm @ 10 psig
Max flow rate	8.0 lpm	9.3 lpm
Thermal overload protector	Thermal fuse; functioning temperature 150°C	
Weight	1.5 kg	
Size (l x w x h)	165 x 165 x 108 mm	
Noise level	58 ±3 dBA	

\*data on file at Philips Respironics

Respironics, InnoSpire, Essence and SideStream are trademarks of Koninklijke Philips N.V. and its affiliates, all rights reserved.



Scan with your smartphone or tablet to go directly to our website.

Please visit [www.philips.com/innospire](http://www.philips.com/innospire) for more information



© 2014 Koninklijke Philips N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

BROUDY HT 5/17/14 MCI 4106223 PN 1103686

Philips Healthcare is part of Royal Philips

Philips Respironics Asia Pacific  
+65 6882 5282  
Philips Respironics Australia  
+61 (2) 9947 0440  
Philips Respironics United Kingdom  
+44 800 1300 845

Respironics Respiratory Drug Delivery (UK) Ltd  
Chichester Business Park, City Fields Way  
Tangmere, Chichester PO20 2FT UK

Customer Service  
+1 724 387 4000  
[www.philips.com/respironics](http://www.philips.com/respironics)

